1<sup>st</sup> EDITION | WORLDWIDE



# THE GLOBAL NAME IN VETERINARY DENTISTRY

DISCOVER THE iM3 DIFFERENCE

WE ARE COMMITTED TO PROVIDING VETERINARIANS WITH THE VERY BEST PRODUCTS AND AFTER SALES SUPPORT WORLDWIDE

VET DENTAL UNITS VET DENTAL X-RAY ULTRASONIC SCALERS HAND INSTRUMENTS

ACCESSORIES & CONSUMABLES

## Table of contents

Hand Instruments	Chapter 1
	Page
Rabbit and Rodent	4-7
Table Top Gag and LED Light	7
Soft Tissue Protector and	
Diamond Disc Safety Shield	8
LED Advantage Straight Nose Cone	9
Canine & Feline Instruments	10-14
Extraction Forceps	15
Feline & Small Dog	16
Curved Deciduous Elevators	17
5pc and 14pc Instrument Sets	18
4pc and 6pc Winged Elevator Sets	19
12pc Extraction Set	20
Feline 8pc Instrument Set	20
Mouth Gags	21
Instrument Tray	21
Safety Glasses	21
Lip Retractor and Oral Curettes	22
Dental Mirror	22
Clear Dental Models	23

#### Consumables

	Page
Dental Burs	24-27
e-VDS Dental Scoring System	28
SANOS Dental Sealant	29
iC Plaque Disclosing Solution	30
CLS Solution	30
Disposable Prophy Angles	31
Prophy Cups	31
Prophy Paste	31

Handpieces	Chapter 3
	Page
High Speed Handpieces	32-33
iPush Tool	33
Low Speed Handpiece & Contra Angles	34
Straight Nose Cones	35
Lubrication Oil	35

Visit our NEW website at www.im3vet.com

iM3 is pleased to Launch its new worldwide website **www.im3vet.com.** Veterinarians now have the opportunity to order direct from iM3 via the website, view the latest product catalogue, subscribe to informative iM3 newsletters featuring case studies and learn more about the iM3 product range.

# Visit the iM3 website and register your email to receive more information.

have been using iM3 products in a veterinary dental referral practice for many many years. They sell top quality products, and have the very best customer service in the industry.

They make the lives of us in the veterinary dental business much easier! We are able to get our jobs done much more efficiently and proficiently.

Kudos iM3! Keep up the great work

Erik Christensen, CVT Arizona Veterinary Dental Specialists





#### Ultrasonic Scalers

#### Chapter 4

	Page
iM3 42-12	37
Ferrite Rod and Scaler Tips	37
iM3 P6 & P6 LED Piezo	38
iM3 P4 & P6 Scaling Tips	38
Cavitron Scaling Tip	38
iM3 Ultra LED	39

#### Dental Units

#### Chapter 5

	Page
Vet-Tome (automated-periotome)	40-41
LED Curing Light	42
iM3 PL3 Polisher (micromotor)	43
iM3 SP6 Scaler & Polisher	43
iM3 GS Deluxe START	44
iM3 GS Deluxe	45
iM3 GS Deluxe LED	46-47
iM3 elite LED	48-49
iM3 Pro 2000 LED	50-51
Features and Benefits	52
Features and Options	53
Pro X Dental Treatment Centre	54-55
Technical Support & Warranty	56

#### X-Ray

Chapter 2

#### Chapter 6

	Page
Dental X-Ray Positioning Kit	57
Bisecting Angle Guide	57
Digital Dental X-Ray - CR7	58
Revolution 4DC X-Ray Generator	59
Hand Held X-Ray Generator	60


#### Practical Products Chapter 7

	Page
iM3 Practical Products	62
About iM3	63



## **Rabbits and Rodents**

iM3's rabbit and rodent instruments have been developed with the help of leading veterinary exotic specialists.

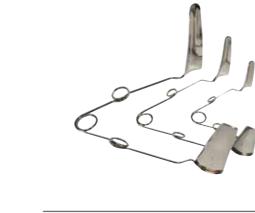

Cheek Di

Cheek Di

#### Cheek Di

The long reach molar teeth. Th the cheeks aw help of Peter I Association (B

## Rodent W Lightweight de rats and small i



Rodent M Designed by i plastic reducir weight on the

**G** The coa I definitely wo beautifully sm Judy Rochette





Cat. No D1026 Used to aid extraction of rabbit incisors. The instrument is used to break down periodontal ligaments on the distal and mesial aspects of the incisor teeth.

D1027

D1028

D1031

Molar Luxator
Designed for the extraction of rabbit cheek teeth. This instrument is used to break down periodontal ligaments on all four aspects of the premolar and molar teeth.

Spatula
Enhances examination of the teeth by repositioning the tongue and cheeks.
Features a light weight handle.

Cheek Dilator – X-Small Suitable for very small mammals where the small cheek dilator is found to be too big.	D1035
Cheek Dilator – Small	D1030

Ch Suitable for smaller mammals such as chinchillas and guinea pigs, enhancing inspection of the oral cavity.

Cheek Dilator – Large Suitable for larger rabbits.

Incisor Luxator











D1046 and D1033 in use.



ilator Long Reach Curved – Small 28mm	<b>Cat. No</b> D1146
ilator Long Reach Curved – Medium 40mm	D1046
ilator Long Reach Curved – Large 60mm	D1047
h 40mm cheek dilator is ideal when using long diamond burs to reduce he radius of the paddle provides better access and viewing as it pushes vay from the molar teeth. iM3's unique design was developed with the Haseler B.Vet.Med. MRCVS. Past president of the British Veterinary Dental VDA) Member of BVDA and EVDS. (not included in D1034 8pc kit)	
<b>Vire Mouth Gag</b> lesign with light spring action. Suitable for chinchillas, guinea pigs, rabbits. Enhances inspection of oral cavity.	D1032
<b>Nouth Gag</b> M3. The Jaws of the new Rodent mouth gag are covered in soft ng trauma. This gag is lighter than comparable gags to reduce the e animals jaws. Autoclavable.	D1033
ated gags are very "kind" to the teeth and soft tissues - ouldn't go back to uncoated gags again. Your gags have a nooth screw mechanism". e DVM, FAVD, Dipl AVDC	
Soft plastic coating reducing trauma	

#### **Rabbit Molar Cutters**

Ideal for trimming excessively long cheek teeth which should then be smoothed off using the extra long diamond bur (Pictured page 8 D2235). By reducing the length with the molar cutters first, this will reduce bur wear and possible thermal injury to the tooth. Care should be taken to avoid tooth fractures and soft tissue damage.

D1045

| 1 - 5



Cheek Teeth Extraction Forceps (Molar Teeth) Specifically designed for the removal of premolar and molar teeth in rabbits. The forcep jaws are narrow and angled at 100 degrees to improve access to the cheek teeth.



Cat. No

D1029



TTMG - Table Top Mouth Gag (Rabbit/Rodent) The TTMG was designed by iM3 with the help of specialists to hold the anaesthetised rabbit or rodent in the best possible position for examination and treatment of the oral cavity without assistance.

- Optional Heater Pad

Table Top Mouth Gag (Rabbit/Rodent)







The TTMG is also suitable for oral surgery!



LED Light

The LED light is ideal for illuminating the oral cavity during rabbit dentistry. The flexible 30cm neck allows the bright light to be easily focused on the area of interest. The LED light is also suitable for the consulting room and surgery where a bright white light is needed. The LED light is supplied with a 110 -240v transformer and 3m long power cord.



D1075 in use - Lower Incisors



D1076 in use - Upper Incisors



Incisor Luxator for labial and lingual/palatal tooth aspects The luxators are used to cut the periodontal ligament and expanding the alveolus by inserting the sharp instrument into the socket. Both incisor luxators are needed to extract the upper and lower incisors.





8 Piece Complete Rabbit and Rodent Kit 8 piece complete set in zipper case . Each instrument is individually available.	<b>Cat. No</b> D1034
Contains:	
Incisor Luxator	D1026
Molar Luxator	D1027
Spatula	D1028
Cheek Dilator - Small	D1030
Cheek Dilator - Large	D1031
Rodent Wire Mouth Gag	D1032
Rodent Mouth Gag	D1033
Cheek Teeth Extraction Forceps (Molar Teeth)	D1029
This instrument set was developed with the help of Dr David Vella BSc BVSc (Hons)	

Diplomate ABVP (Exotic Companion Mammals).

| 1 - 6

Hand Instruments

## iM3 Small Animal Dental Equipment

- Adjustable height stainless steel gags to suit any size rabbit or rodent
- Adjustable table height ensures the best viewing angle of the oral cavity
- Washable silicon pad under the patient reduces heat loss and is easy to clean
- The complete TTMG is easy to clean and autoclavable
- Quality components and innovative design for ease of use.

• Temperature controlled silicon heater pad further minimises body temperature loss. Supplied complete with low voltage power supply with international approvals. • LED Light for improved illumination of the oral cavity. See T2200 below.

Orange heated silicon pad & transformer



1 - 7



1 watt bright white LED light with 30cm flexible neck.



Complete set

Soft Tissue Protector and Diamond Disc Safety Shield

Cat. No

Dental Instruments for rabbits, guinea pigs and chinchillas. The iM3 soft tissue kit is an ideal set of instruments for reducing and leveling molar and incisor teeth.


NEW

Complete set	D2245
<ol> <li>Diamond Disc Cutter</li> <li>22 mm, coated on both sides with HP mandrel</li> </ol>	D2200
2) Safety Shield Fits L6860. Safety shield for diamond disc, stainless steel	D2220
<ol> <li>Soft Tissue Protector Fits L6860. The soft tissue protector guards the tongue and cheeks by exposing only the bottom cutting surface of the diamond bur when used on cheek teeth.</li> </ol>	D2230
<ul> <li>4) Long Cutting HP Bur Diamond bur X-Long 60 mm. Diamond cleaning stone available (D9750)</li> </ul>	D2235
The above items fit a straight nose cone 1:1 (Blue) fitted to a low speed motor. (shown below).	
<ol> <li>Straight Nose Cone 1:1 (Blue)</li> <li>Straight Nose Cone 1:1 ratio (blue) use HP burs for rabbit and orthopaedic work. Autoclavable.</li> </ol>	L6860
All parts are available separately. 5	Constant Street



Light in the mouth!

We purchased your LED Advantage Straight Nose Cone and the L6200 Low Speed Motor to fit our countertop dental unit. This nose cone has increased our ability to perform dental procedures on our smaller pocket pets. The LED light allows the veterinarian to easily visualize the molar teeth and the varying sizes of burs help to rapidly improve the occlusion on the cheek teeth. It has dramatically cut the time needed for these procedures. We couldn't be happier with this unit and recommend it to anyone doing pocket pet dentistry!" Mary L. Berg BS, RLATG, RVT, VTS(Dentistry)

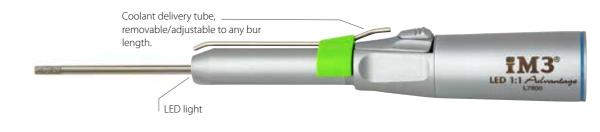

#### Leporidae Skull Model

iM3's new Leporidae skull model allows you to observe the tooth and root pattern of a typical rabbit and demonstrate to your clients why a procedure is necessary and what will be involved. Also included with the model is an education card. The card is double sided and shows the location of the nasolacrimal duct. The other side shows teeth location, name and numbering for easy identification.



Model size 9cm x 4.5cm x 5.5cm (3 1/2" x 1 3/4" x 2.2") Card size 16.5cm x 13cm (6 1/2" x 5 1/4")



HP Burs suitable for rabbit dentistry Fits any straight nose cone HP fitting. Long cutting HP Diamond Bur Diamond bur X-long 60mm, can be used with or without D2230 soft D2235 tissue protector. The diamond bur is less likely to cut soft tissue. Narrow Conical HP Bur for Rabbits The carbide bur is suitable for levelling molar teeth. Care is needed with D9550 the soft tissue. Does not fit D2230. Soft Tissue Protector. Large Tubular HP Bur for Rabbits The large carbide bur is suitable for levelling molar teeth. Care is needed D9555 with the soft tissue. Does not fit D2230.

#### iM3 Small Animal Dental Equipment

#### 1:1 LED Advantage Straight Nose Cone

The 1:1 LED Advantage fits any low speed motor or electric micromotor • Self generated LED white daylight, produced by the drive shaft of the low speed air motor or electric micromotor

Large illumination area providing outstanding visibility

- Coolant delivery tube, removable/adjustable sterile coolant to any bur length
- Ideal for rabbit dentistry and orthopaedic work. Autoclavable.
- Takes HP burs and 3/32" drills

Cat. No D1070

PORIDAE SKULL time in the loss her IM3

L7800

| 1 - 9

Hand Instruments

# iM3 Luxators and Winged Elevators

The dental instruments that meet the needs of the modern day Veterinarian

iM3 Luxators are very fine and sharp. The 5 sizes make them ideal for both canine and feline teeth. They are used and recommended by dental specialist's for cutting the periodontal ligament and expanding the alveolus (tooth socket).

iM3 winged elevators are available in 6 sizes and are used to elevate and extract the tooth. The winged shape of the elevator tip matches the curve of the tooth and the extra leverage aids in extraction.

vetering specialists are developed with the guidance of interiments are developed with the guidance of a lifetime waranty\*







Periosteal E Used for elevating Periosteal Elevato



\*Lifetime warranty covers faulty material and workmanship. Damage due to the incorrect use or handling of the instrument is not covered by the lifetime warranty.


Veterinary Winged Dental Elevators	
iM3 winged elevators are light-weight yet extremely durable. The sharp, tapered winged tip allows easy insertion into the PDL. The width of the winged tips have been shaped to match the size and shape of the tooth's root structure (sizes 1-6 mm). This unique winged shape increases control of the instrument and provides additional leverage to aid in extraction of the tooth.	
d Elevators with Standard Handle	Cat. No
evator 1mm	D1001
evator 2mm	D1002
evator 3mm	D1003
evator 4mm	D1004
evator 5mm	D1041
evator 6mm	D1042
oot Tip Pick	
g with extraction of retained root tips or root fragments. pick - stubby handle	D1005 D1005S
Elevator Ig the muco-periosteal flap during extractions on large breed dogs.	D1006
or - stubby handle	D1006S

11-11

	iM3 Winged Elevators with Stubby Handle The stubby handle elevators available as a 6 pc. set. See page 16	Cat. No
	Winged Elevator 1mm	D1001SH
Harrison and Dispatch	Winged Elevator 2mm	D1002SH
	Winged Elevator 3mm	D1003SH
des industries	Winged Elevator 4mm	D1004SH
	Winged Elevator 5mm	D1005SH
	Winged Elevator 6mm	D1006SH

6 pc Stubby Handle Set includes stainless steel instrument holder/case

D1044

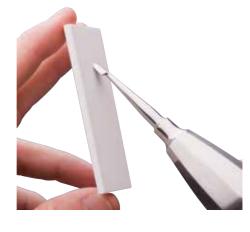

iM3 Wedge Dental Sharpening Stone The wedge sharpening stone has a flat surface with a radius on either side. Dental hand instruments should be kept sharp and true to their original design. For best results, instruments should be sharpened after use. This Arkansas stone should be used with water or oil. Sharpening instructions included.

A3110







Dr Loic Legendre at the AVDS dental workshop. Sponsored by iM3

1 - 12

iM3 Luxators, sizes 1-5mm are very fine and sharp. The 5 sizes allow their use for both feline and Canine teeth. They are used and recommended by dental specialist's for cutting the periodontal ligament and expanding the alveolus (tooth Socket) and are ideal for

The stubby (shorter length and larger diameter) handle makes extraction so much safer as the users fingers can be placed at the instrument tip. Luxators are not elevators, we recommend winged elevators are used to elevate and extract the tooth.

#### 5pc Luxator instrument set (1-5mm Stubby)

D4501 1mm Luxator (Stubby Handle)

D4502 2mm Luxator (Stubby Handle)


D4503 3mm Luxator (Stubby Handle)

D4504 4mm Luxator (Stubby Handle)

D4505 5mm Luxator (Stubby Handle)

Cat. No D4500 m

È

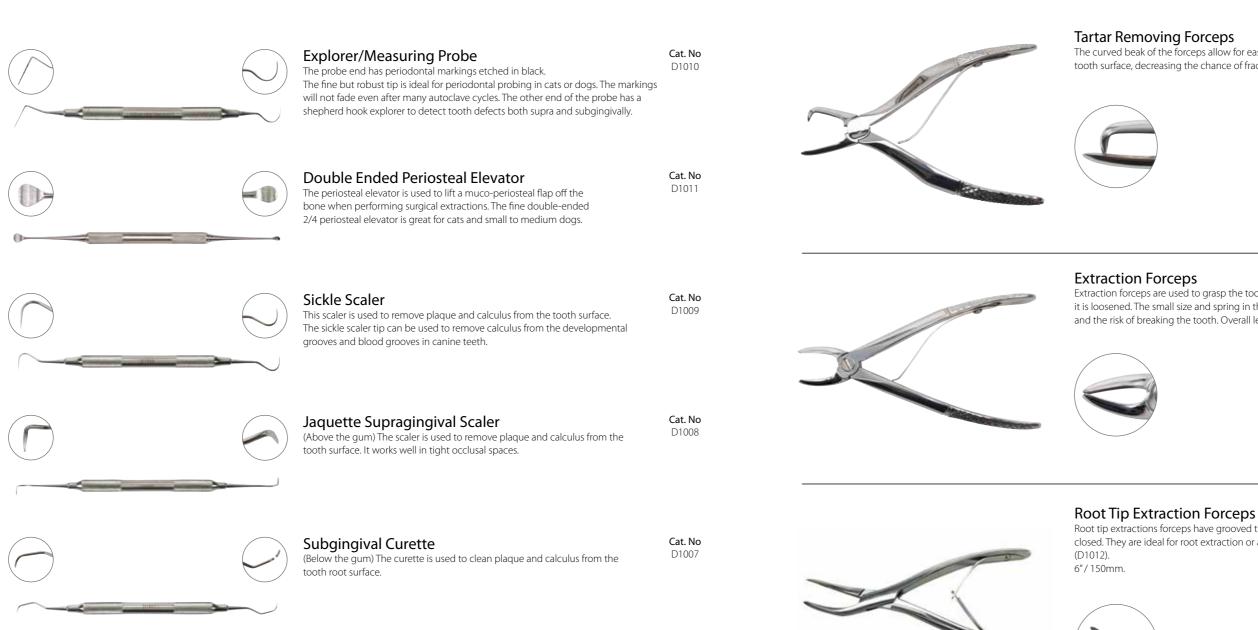



#### Minnesota Retractor

Used to lift and protect the tissue and the user when performing surgical extractions whilst operating the high speed drill. One end holds the cheek back and the other holds the tongue down, a multi purpose retractor.



Cat. No D1014



5 1/2"/ 138mm.



| 1 - 14

Hand Instruments

The curved beak of the forceps allow for easy removal of heavy calculus from the tooth surface, decreasing the chance of fracturing the tooth.

Cat. No D1013

Extraction forceps are used to grasp the tooth and remove from the socket after it is loosened. The small size and spring in the handle helps prevent excessive force and the risk of breaking the tooth. Overall length 125mm (5").

Root tip extractions forceps have grooved tips with fine curved ends that meet when closed. They are ideal for root extraction or as an alternative to regular extraction forceps

Cat. No D3050

Fragment Extraction Forceps

Root fragment forceps are ideal for extracting mobile broken root fragments and small, delicate teeth. The grooved and serrated jaws assist with grasping fragments. Cat. No D3060

Cat. No D1012 1 - 15



Instruments mentioned on this page are available as a set (Cat. No D2000)



Cat. No

D2001

Cat. No

D2002

Cat. No

D2003

Cat. No

D2005

Cat. No D2006

Cat. No

D2007


Cat. No

D2008

Cat. No

D2004

encountered.

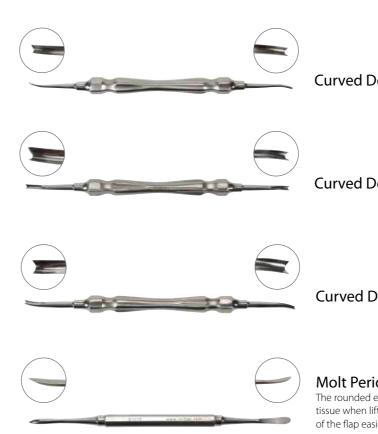



| 1 - 16

#### Curved Deciduous Elevators

The curved deciduous elevators are ideal for the removal of retained canines (milk teeth) and are available in three sizes suitable for the different sized deciduous teeth

The sharp elevator end cuts the PDL whilst the curve of the elevator tip follows the mesial curve and the other end the distal curve of the root. Curved deciduous elevators make the removal of retained teeth quick and easy because they follow the curve of the root.

eciduous Elevator – Small (2mm)	<b>Cat. No</b> D1016
eciduous Elevator – Medium (3mm)	Cat. No D1017
eciduous Elevator – Large (4mm)	<b>Cat. No</b> D1018
osteal Elevator – No. 9 dge of the periosteal instrument minimise the tearing of ting the gum flap. The sharp end makes the initial lifting er. Suitable for medium to large dogs.	<b>Cat. No</b> D1015

## Small Dog or Cat Right Angle Extraction Forceps

Designed by iM3 for easier extraction of cat, small dog and deciduous teeth. Parallel right angle jaws reduce the risk of tooth fracture during extraction. By not meeting at the end, they reduce the risk of breaking off cats teeth. Small handles provide a more tactile feel.

Cat. No D1019

1 - 17

Instrument Sets

iM3's Instrument sets take the guess work out of choosing the correct dental instrument. iM3 instruments are used throughout the world in teaching

iM3 5 Piece Specialised Instrument Set The set is intended as a supplement to ordinary tooth extraction sets.

The set consists of: Molt Periosteal Elevator - No 9 Curved Deciduous Elevator - small Curved Deciduous Elevator - large

#### Cat. No

iM3 14 Piece Instrument Set, Stubby Handle Elevators Cat. No D1100

The set consists of: Winged elevators Size 1-5mm Straight root tip pick Tartar Removing Forceps Extraction Forceps Minnesota Retractor Subgingival Curette Jaquette Supragingival Scaler Sickle Scaler Explorer / Measuring Probe Double ended Periosteal Elevator Each instrument is available individually

# The curved deciduous elevators are ideal for the removal of retained canines (milk teeth). Each instrument is available individually.

Curved Deciduous Elevator - medium Extraction Forceps, right angled jaw

#### iM3 14 Piece Instrument Set, Standard Handle Elevators D1000

The set consists of: Winged elevators Size 1-5mm Straight root tip pick Tartar Removing Forceps Extraction Forceps Minnesota Retractor Subgingival Curette Jaquette Supragingival Scaler Sickle Scaler Explorer / Measuring Probe Double ended Periosteal Elevator Each instrument is available individually



Stubby Handle Set



## iM3 6 Piece Winged Elevator Set Each elevator is available individually.

Standard Handle Set



#### iM3 6 Piece Winged Elevator Set

Stubby Handle Set



## iM3 Small Animal Dental Equipment





#### iM3 4 Piece Winged Elevator Set

iM3 Winged Elevator sizes 1-4 mm in standard handle, supplied in a black pouch. Each elevator is available individually.

Cat. No D1038

Standard Handle Set

iM3 4 Piece Winged Elevator Set

iM3 Winged Elevator sizes 1-4 mm in stubby handle, supplied in a black pouch. Each elevator is available individually.

Cat. No D1039

iM3 Winged Elevator sizes 1-6 mm in standard handle, supplied in a black pouch.

## Cat. No

D1048

iM3 Winged Elevator sizes 1-6 mm stubby (shorter length and larger diameter) handle. The stubby handle makes extraction so much safer as the user's finger can be placed at the instrument tip. The new handle design is ideal for small or large hands. Supplied in a black pouch. Each elevator is available individually.

Cat. No D1049

6

Ļ



#### iM3 6 Piece Winged Elevator Set

iM3 Winged Elevator sizes 1-6 mm stubby (shorter length and larger diameter) handle. The stubby handle makes extraction so much safer as the user's finger can be placed at the instrument tip. The handle design is ideal for small or large hands. The set is supplied in a stainless steel autoclavable instrument holder/case. Each elevator is available individually.

Cat. No D1044

Stubby Handle Set incl. stainless steel instrument holder/case

#### iM3 12 piece Extraction Set - Updated

The complete extraction set provides the veterinarian with the instruments required to perform both simple and surgical extractions.

#### The set consists of:

D1001SH - D1006SH (1-6mm Stubby handle winged elevators) D1011 Periosteal elevator double ended 2/4 D1012 Extraction Forceps Universal T1000 Needle holders – Mayo Hegar T1200 Tissue forceps – Adson Brown T1300 La Grange Scissors, curved, one serrated blade T1400 Scalpel blade holder with round ergonomic handle D1043 Supplied in stainless steel autoclavable instrument cassette. 20.5 x 12 x 3cm

## NEW

iM3 13 piece Extraction Set


The set consists of the same 12 instruments listed above and in addition to : T1500 Curved Blunt Suture Scissors The 13pc set is supplied in a large stainless steel autoclavable cassette 20.5 x 16.5 x 3cm

#### iM3 Feline 8 Piece Instrument Set

The feline periodontal kit has been designed exclusively for use in cats and small dogs. Its design allows the practitioner to perform a complete oral examination (including detection of feline resorptive lesions) and cleaning of all tooth surfaces - both supra and subgingivally. All instruments are available individually.

Cat. No D2000

The set consists of: 1 mm Wide Straight Tip Elevator, 2 mm Wide Straight Tip Elevator Probe/Explorer Extraction Forceps Periosteal Elevator Resorptive Lesion Probe Sickle Scaler and Curette.

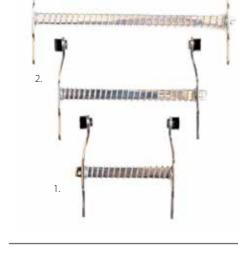

# iM3 Mouth Gag - 12 cm opening (3) iM3 Mouth Gag 3 piece set












Anti-Fog Safe Your sight is one of protection. Anti foo veterinary dental pr fragments.

www.im3vet.com

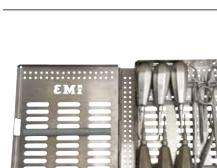

ents

Hand Instrum

D4100 Pictured

	HM3	1
-	-	
	-	
	-	
	1.10	-

Cat. No D4000

> Cat. No D4100

1



#### iM3 Stainless Steel Mouth Gags

Unlike other (heavy duty spring-loaded) mouth guards that tend to overextend the temporomandibular joint, the iM3 designed mouth guard has a unique locking system that requires little force to hold the gag at the desired open position inside the animal's mouth. The lightweight gag and spring coupled with the unique locking system reduces the risk of myalgia, neuralgia, and/or trauma to the temporomandibular joint.

iM3 Mouth Gag - 5 cm opening (1)

iM3 Mouth Gag - 8 cm opening (2)

Don't confuse our Veterinary specific mouth gags for these heavy duty monsters!

Needle Holders with Tungsten Carbide Jaws
Cardide Jaws
Designed for multiple tasks and ideal for wiring jaws or
orthopaedic work.
The strong but fine jaws have tungsten carbide inserts
for wire cutting.

#### Scalpel Blade Holder - Round 13cm (5<sup>1</sup>/<sub>4</sub>")

The iM3 scalpel holder with tapered round handle is ergonomically designed (shaped like a pen) providing greater control during use. Size 3 handle.

Autoclavable Instrument Trays Autoclavable to 135° C. Instrument Tray 95 mm x 195 mm (1)	<b>Cat. No</b> D1020
Instrument Tray 340 mm x 245 mm. Has 10 drilled bur holes (2)	D1036
Stainless steel instrument holder/case 205 mm x 120mm x 30 mm (Excl. instruments) (3)	D1043
Large Stainless steel instrument holder/cassette 205 x 165 x 30mm (4)	D1096

fety Glasses	Cat. No
f your most precious assets and deserves the best possible	G0000
g safety glasses will help protect your eyes during routine	
procedures from aerosols as well as tooth or calculus	

Cat. No

D1023

D1024

D1025

D1040

Cat. No

D3000

Cat. No

T1400



Feline Clear Jaw Model Card Size 16.5 cm x 13 cm (6 1/2" x 5 1/4")

will be involved.





<b>Oral curettes</b> D1086 Oral Curette double ended small (Feline & Canine)	<b>Cat. No</b> D1086
D1087 Oral Curette double ended Med - Large (Canine)	D1087
The two sizes of Oral Curettes (spoon shaped) are ideal for debrid- ing the alveolus socket after tooth extraction.	

Cat. No D1095

For large jowled breeds, using 2 lip retractors makes seeing your work site a breeze!



DoxiProbe The new DoxiProbe was designed (by Drs. Paul Mitchell and Randi Brannan, both AVDC diplomates) to aid in detecting and treating periodontal pockets. The probe measures pocket depth and the lattened, paddle-like end helps place Doxirobe in the prepared pocket.



Canine Clear Jaw Model Card size 16.5 cm x 13 cm (6 1/2" x 5 1/4")





www.im3vet.com

1 - 22



#### Leporidae Skull Model

iM3's new Leporidae skull model allows you to observe the tooth and root pattern of a typical rabbit and demonstrate to your clients why a procedure is necessary and what

Also included with the model is an education card. The card is double sided and shows the location of the nasolacrimal duct. The other side shows teeth location, name and numbering for easy identification.



Model size 7 cm x 5 cm x 3 cm (2 3/4" x 2" x 1 1/8") Included with the model is a plastic educational card showing the tooth and root location as well as the triadan numbering system for easy identification.



## Model size 12 cm x 7 cm x 5 cm (4 3/4" x 2 3/4" x 2)" Included with the model is a plastic educational card showing the tooth and root location as well as the triadan numbering system for easy identification.



Set includes both the Feline and Canine models pictured above.

Cat. No


D1060

Cat. No D1070

| 1 - 23

Cat. No D1050

Cat. No D1055

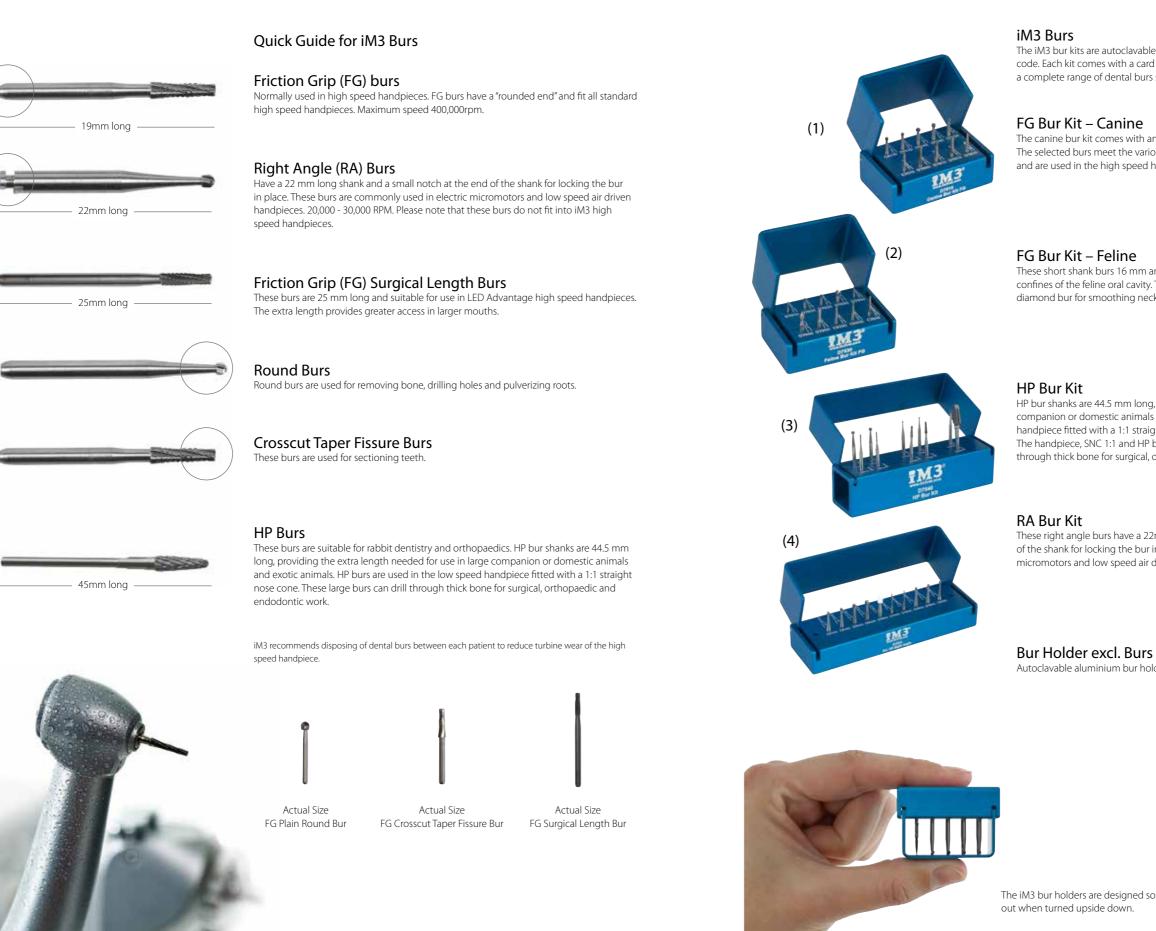

www.im3vet.com

<b>iM3 Burs</b> The iM3 bur kits are autoclavable and have an illustration of the bur head and reorder code. Each kit comes with a card suggesting the use for various burs in the kit. iM3 offers a complete range of dental burs sold individually in packets of 5 burs.	
FG Bur Kit – Canine The canine bur kit comes with an assortment of FG burs 19mm long. The selected burs meet the various requirements of a veterinary clinic and are used in the high speed handpiece. (1)	<b>Cat. No</b> D7510
FG Bur Kit – Feline These short shank burs 16 mm are preferred for use in the confines of the feline oral cavity. This kit contains an extra fine diamond bur for smoothing neck lesions (CLL's, FORL's). (2)	D7530
HP Bur Kit HP bur shanks are 44.5 mm long, providing the extra length needed for use in large companion or domestic animals and exotic animals. HP burs are used in the low speed handpiece fitted with a 1:1 straight nose cone The handpiece, SNC 1:1 and HP bur can be autoclaved. These large burs can drill through thick bone for surgical, orthopaedic and endodontic work. (3)	D7540
<b>RA Bur Kit</b> These right angle burs have a 22mm long shank and a small notch at the end of the shank for locking the bur in place. These burs are commonly used in electric micromotors and low speed air driven handpieces. 20,000 - 30,000 RPM. (4)	D7520

Autoclavable aluminium bur holder for 10 FG burs. (5)

(5)

The iM3 bur holders are designed so that the burs will not fall


D7500

			Bur Type		Bur Type Head Size (		Cutting	Friction Grip for High Speed Handpieces			SW Handpieces
Dental Bur Chart			гтуре	mm)	Length (mm)	FG Feline (16mm)	FG Standard (19mm)	FG Surgical (25mm)	RA (19mm)		
	high speed handpieces. iM3 recommends disposing										
s between each patient to reduce tu	rbine wear of the high speed handpiece.		1/4	005	-	D7910	D7911	D7915	C9091		
·			1/2	006 008	-	D7909	D7912 C9050	D7916	C9092 C9093		
1 1 9 0			2	010	-	D7908	C9080	C9085	C9095		
	FG Plain Round Bur	Plain Round	3	012	-	27,500	0,000	0,000	C9096		
	Pack of 5 autoclavable burs.		4	014	-	D7920	C9110	C9115	C9100		
	For removing bone, drilling holes and pulverizing roots.		5	016	-				C9097		
	iM3 recommends size 2 & 4 for cats and small dogs and		6	018	-	D7907	C9140	C9145	C9160		
ACTUAL SIZE	size 4 & 6 for larger dogs.		8	023	-	D7930	C9142	D7935	C9094		
1/4 1 <b>2 4 6</b> 8			699	009	4.2		D8012				
			700	010	4.2	C9001	D8013	D8015	C9098		
7910 C9050 C9080 C9110 C9140 D7930			701	012	4.2	C9002	C9320	D8016	C9330		
		Cross-cut Fissur Taper	701L 702	012 016	6.0 4.4		C7010 C9350	C7015	C9390		
			702 702L	016	6.0		D7810	D7815	CA3A0		
			702L	010	4.8		C7030	01010			
	FG Crosscut Taper Fissure Bur		703L	021	7.5		C4450				
	Pack of 5 autoclavable burs.		556	009	4.2		C9279	C9275			
	For sectioning teeth.		557	010	4.2		C9280		C9290		
	iM3 recommends size 699 & 701L for cats and small dogs and	Cross-cut Fissur		010	6.0		C9278				
	size 702L for larger dogs.		558	012	4.2		C9250	C9255			
ACTUAL SIZE			559	014	4.4		C9300				
<b>99</b> 701 <b>701L</b> 702 <b>702L</b> 703 703L			169	009	4.2		D8011	D0014			
		Fissure Taper	169L 170L	009 010	5.2 6.0		D8006 D8009	D8014 C9089			
3012 C9320 C7010 C9350 D7810 C7030 C4450			170L 171L	010	6.0		D8009	C9009			
			171L	012	6.0		D8003				
		Round-end	1701	012	4.2		D7949				
		Cross-cut Fissur Taper	1702	016	4.4		D7950				
	FG Surgical Length Burs	Round-end	1557	010	4.2			D7801			
	Pack of 5 autoclavable burs.	Cross-cut Fissur	1550	012	4.2			D7802			
	The extra length provides greater access in large mouths.	Round-end Fissu		012	4.2		D7800				
	These burs are 25mm long and only suitable for use in the		33 1/2	006	0.5	D7940	C9229	D7803			
	Advantage high speed handpiece.	Inverted Cone	35	010	1.1		C9230		07020		
			37 	014 006	1.4 1.2		C9231 D8329		D7830		
ACTUAL SIZE			330	008	1.6		C9440				
1/2 2 4 6 701L 702L		Pear	331	010	1.7		D8331				
			332L	012	4.2		D8332				
7916 C9085 C9115 C9145 C7015 D7815			333L	014	4.4		D8333				
			Assorted	-	-		D7560				
			Round 4	014	-		D7040				
			Round 6	018	-		D7060				
	iM3 FG Diamond Assorted Bur set – 5pc	Cat. No	Round 8	023	-		D7080				
A DE LE		D7560 Diamond	Taper Flat End Cylinder	014 016	8.0 8.0		D7091 C9400				
9 9 0			Flame 861GK	016	6.00		D7820				
	iM3 medium grit (blue band) diamond bur set featuring precision European		Football (016)	012	3.5		D7960				
	made diamond burs. Suitable for use in any high speed dental handpiece.		Football (023)	023	5.0		D7095				
	The set consists of 1 x size 4, 6 and 8 plain round and 1 x round end taper and		Egg 379	023	4.2		D8379				
	football bur.		Arkansas Flame #661	025	7.0		D7851		D7850		
		Others	Metal Cutter - 557E	010	4.2		D8557				
	All iM3 diamond burs are available individually in a pack of 5. See order codes below.		Flame 7903	012	3.6		D8010				
ACTUAL CITE	שבי טועפו נטעפג שפוטא.	The	100	78	W MP	0		₩.			
ACTUAL SIZE		14		1				11	14		
4 6 8 Taper Football			1.					1.7	11		
07040 D7060 D7080 D7091 D7095		17	113	11				17	11		
		19	100					100	100		

2 - 26

# Size **699** 701 **701L** 702



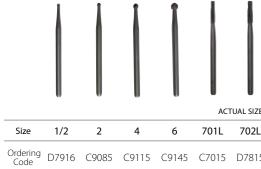



Size	4	6	8	Taper	Football
Ordering Code	D7040	D7060	D7080	D7091	D7095

U
نب
Ū
>
$\mathbf{m}$
Ξ
• —
>
>
$\leq$
$\geq$

## iM3 Small Animal Dental Equipment

Introducing

# HM3

**E-VDS**Plus electronic Veterinary Dental Scoring by the European Veterinary Dental Society (EVDS) The revolutionary dental scoring system that's just a click away!

Proudly sponsored by iM3 The Veterinary Dental Company & Virbac.

## How does it work?

- Ist step: registration as a user
- 2nd step: setting up clinic data
- 3rd step: selecting language for the printout; five standard languages and an option for translation into any language
- Scoring on dental chart by clicking the mouse
- Choice between puppy or adult dentition and single or multi tooth scoring
- Recorded data can be printed, downloaded as PDF or sent by email

## The value of e-VDS Plus For the vet:

- The clinical findings are scored on dental charts (by David Crossley) which provides thorough documentation and history for each patient
- The tutorials made available during the process will help train vets in dentistry with minimal effort
- Individual clinic scoring sheet with your own practice logo and header will help build reputation and loyalty

#### For the pet owner:

- The owner is involved in the process
- Printout with the assessment of their pet
- A legend explaining the findings of the scoring
- Provides a chance to follow the development of treatment or prophylaxis
- Encourages follow-up



## Why e-VDS Plus?

iM3 and Virbac's role, as worldwide leaders in veterinary dentistry, is to help the market expand. e-VDS Plus was designed to make dentistry more systematic and encourage vets (or their staff) to spend enough time looking at the oral cavity during the annual check.e-VDS is an interactive way to raise the level of education of general practitioners in veterinary dentistry, thanks to a practical process and unique iconography.

## Where to find it?

You can find e-VDS Plus on the brand new EVDS website www.evds.org or on the iM3 or Virbac website.



The iM3 Dental Unit is superior in every way. The handpieces are of much higher quality, the CLS flush system gives me piece of mind...NO COMPARISON! PRO2000 and 42-12, Carbondale, Co USA



#### iC Plaque

The iM3 iC Plaque Disclosing Swabs show when plaque is present on the teeth. These applicators are quick and easy to use. By simply rolling the swab gently across the teeth at the gum edge the plaque stains pink to indicate germ-filled layers on the teeth. This is a visible sign of periodontal disease. And it's there, simple and easy for your client to see.



## Guaranteed not to catch hair or whiskers! Also reduces heat/friction





# be single use only.

iM3 iClean Prophy Paste

holder supplied.

iClean R&R Prophy Paste, Med Grit, no Fluoride – Bag of 100 x 2g individual units complete with finger holder

## iM3 iClean Prophy Paste, 255 g Tub

A fine grit prophy paste with Fluoride and fruit flavour (smells like bubblegum) iClean R&R polishing paste has been specially formulated for use in animals. The fine grit paste is suitable for all polishing cups. Adequate dental polishing is a must!

iClean Prophy Paste with Fluoride, Fine Grit, Fruit Flavour (smells like bubblegum) – 255g 9oz Tub



#### iM3 iC Plaque, 72/pk

- Plague Indicating Swabs
- Plaque is the cause of periodontal disease
- Plague is invisible until stained
- Highlights the areas where plaque is present you will be suprised by how much plaque is present!
- It is a helpful tool in painting a picture of the animal's oral health for the owner
- been removed during the procedure
- Saves time scaling by indicating plaque and calculus on the teeth
- Will increase dental procedures within the clinic

- Convenient dispenser box



#### CLS Solution 250ml concentrate

CLS is a uniquely formulated dental solution designed to reduce bur cutting time and mask mouth odours (spearmint). CLS contains chlorhexidine and can be used in the water supply of any dental machine or ultrasonic scaler. CLS is guaranteed when used at the correct dilution with distilled water not to block the water lines or coolant system. Easy to mix - just add distilled water to make 1L of working solution.

Important Note:

Products containing chlorhexidine may cause oedema with prolonged oral use. Pack the back of the throat to minimise this risk.

CLS Solution 250ml concentrate - single C3150

CLS Solution 250ml concentrate - 6 pack



- Useful when scaling the animal's teeth to ensure that all the plaque has

- · Quick and easy to use
- Excellent training tool









#### Cat. No H1001

#### iM3 Oscillating Disposable Prophy Angle, 100/pk

Rock and Roll prophy angles are complete, ready to use with a soft rubber prophy cups. The R&R's are pre-lubricated and iM3's unique blue colour is your guarantee of the prophy angle lasting a complete procedure. The 90° reciprocating 'back and forth' action of the prophy cup prevents the animals muzzle hair from being tangled and allows the cup to have continuous contact and pressure on the tooth surface. The reciprocating action also reduces heat and friction on the tooth. Traditional prophy angles spin 360 deg, wind muzzle hair and generate excessive heat.

## iM3 Twister Disposable Prophy Angle 100/pk

The new iM3 twister disposable prophy angles are supplied with a soft ribbed and webbed latex free cup. These premium quality individually packed prophy angles are designed with unique gearing for ultimate smoothness. Guaranteed not to dislodge the cup during use, the cup rotates 360 deg during use. Single use only.

#### iM3 Twister Disposable Prophy cups – Snap on 100/pk.

iM3 replacement disposable snap on prophy cups, can be used on the iM3 L7465 Oscillating Disposable prophy angles. The soft grey cups are latex free and designed to Cat. No P4640

Cat. No L7510

Cat. No 17500

Cat. No L7460

Cat. No

L7465

www.im3vet.com

iClean R&R polishing paste has been specially formulated for use with iM3 R&R oscillating disposable prophy angles. The oscillating action requires a slightly coarser paste than conventional polishing heads. Convenient to use, individual cups help reduce cross contamination and minimise the transfer of viruses between animals. Won't splatter and does not contain flavour oils or fluoride. Packaged in a convenient bag of 100 individual units with finger

2 - 31

Consumables




#### Why iM3 LED light?

With a colour temperature of 5,500 K the iM3 LED Advantage hand piece produces a radiation spectrum that corresponds to daylight. The high colour rendering index (CRI 90) of the LED enables relaxed work thanks to daylight quality light and provides exceptional reproduction of natural colours and clear colour contrasts in the mouth.

## Reliable – Efficient – Precise – iM3 Advantage

Excellent ergonomics

Fatigue -free working thanks to easy 360 degree rotation of the turbine on the coupling, low overall weight and slender design. Unique handpiece design allows component replacement on site by the user

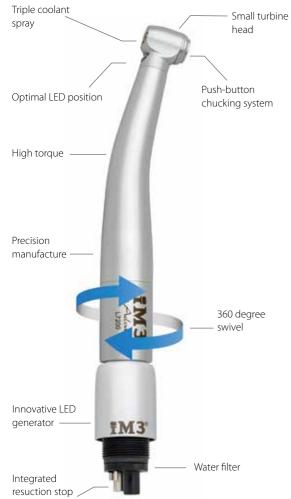

#### LED Advantage Swivel High Speed Handpiece New Midi head provides greater access in small mouths, with all the

Cat. No L7300

- · Self generating LED illumination, no electricity just needs air
- Daylight quality bright white light, unlike yellow F/O
- · Largest illumination area of any dental handpiece

advantages of the iM3 LED handpieces

- Patented turbine design eliminates bacterial/viral suck back
- Autoclavable
- Triple spray coolant outlet for effective bur cooling
- 1 year complete handpiece warranty including the turbine and chuck mechanism
- · Swivel handpiece will reduce fatigue and removes handpiece tension on the operators wrist
- 360° swivel with guick release coupling housing the LED generator
- · Fits all air driven dental machines with four hole connection
- The L7300 includes the L2700 handpiece and the L7400 swivel coupling

"... it seems like such a simple concept, but one that significantly reduces the need to torque your wrist or constantly re-position the handpiece to achieve the appropriate angle of approach to a tooth. It also helps reduce the pull of the weight of the cord on your wrist..."

Lina N. Mohamed, DVM Interbay Veterinary Hospital, Seattle, WA



Easy operation of the push button or iM3 yellow iPush



Precision cooling The triple spray cools the treatment site from all sides




Patented hygenic head system






Replaceable water filter Prevents dirt particles from entering the turbine



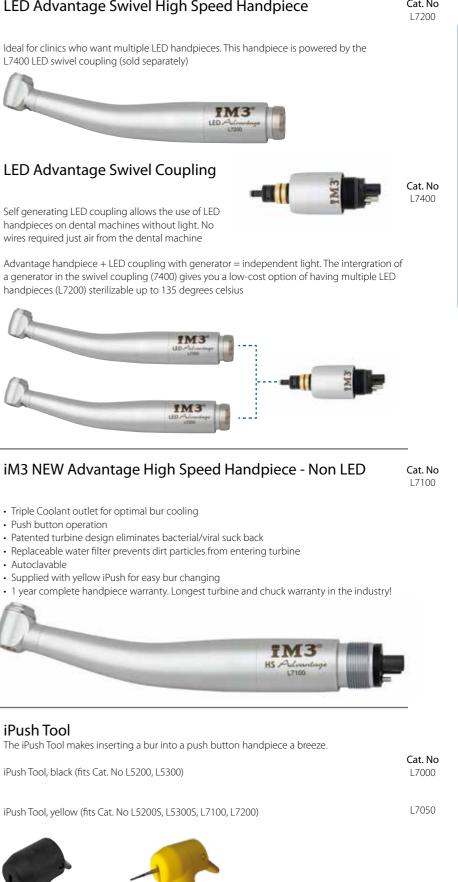

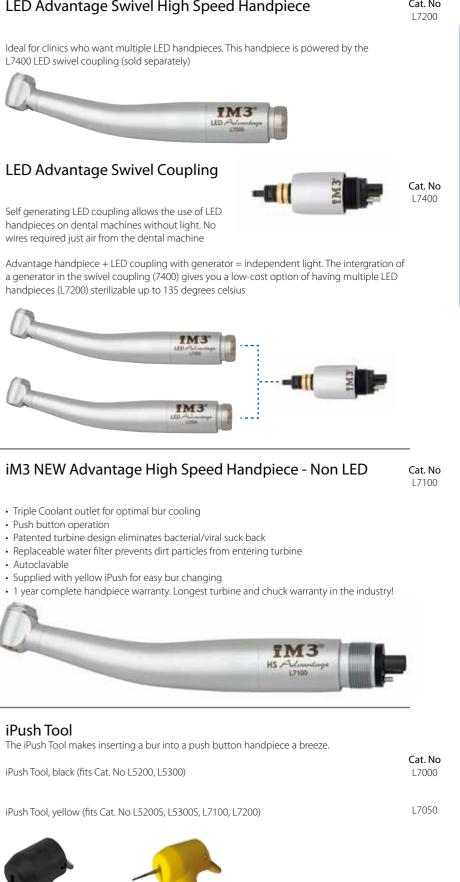










For L5200 and L5300

For L5200s, L5300s, L7100 and L7200 Designed and Manufactured by iM3. EC Patent No. 001882788 USA Patent No 29/398937

Advantage I The motor is Maximum s • Complete	low speed motor er is light weight and s speed 25,000rpm. e 1 year warranty, inc	5	
<ul> <li>Autoclava</li> </ul>		le connection ole forward and reverse also available	L6980
<ul> <li>Reduces ti</li> <li>Does not i is as simpl</li> <li>The Advar</li> <li>Autoclava</li> </ul>	the polishing speed require a chuck ring le as pushing it on!	y. Fitting a R&R oscillating prophy cup t any Low Speed Motor	<b>Cat. No</b> L6250
•	PUSH IN	¢	L6200 + L6250
ales	M3" Lainendage Lazzo	IM3° LS Aufundange	

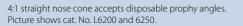

PULL OUT

IM3

Fasy

Friction

Grip









Standard 4-hole connection Internal water supply when connected to a dental machine with LS water. Suitable for running contra-angle and straight handpieces with internal water. Forward/reverse and speed adjustable. Autoclavable 135 Deg C.
4:1 LED Slow Speed Contra-angle
Self generating LED light with internal water. Connects to any air driven or electric Low Speed motor. Single coolant spray, head size diameter 10.1mm 4:1 speed reduction (Green). Push Button chuck system, accepts all RA Burs Use with (L2010) LS motor with internal water. Autoclavable 135 Deg C.
1:1 LED Slow Speed Contra-angle
Self generating LED light and internal water. Connects to any air driven or electric Low speed motor. Single coolant spray, head size diameter 10.1mm 1:1 - no speed reduction (Blue). Push Button chuck system, accepts all RA Burs Use with (L2010) LS motor with internal water. Autoclavable 135 Deg C.
1:1 Slow Speed Contra-angle

Low Speed Motor With Internal Water

360° Swivel

Cat. No

L2010


Cat. No L2002

Cat. No L2004

Cat. No

L2007

o internal water or LED.
onnects to any air driven or electric Low speed motor.
1 - no speed reduction (Blue). Push button chuck system, accepts all RA Burs.
se with standard iM3 LS motor L6200.

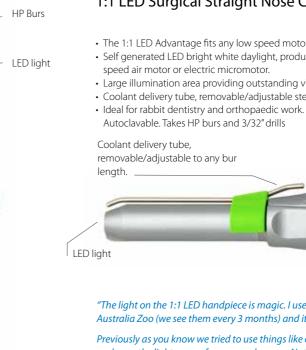



IM3

MD-30

# iM3 Straight Nose Cone 1:1

orthopaedic work. 1:1 means no speed reduction. Autoclavable

## iM3 Advantage Handpie

 Intelligent oil and short • MD-30 is an all syntheti workload of veterinary MD-30 contains no toxid and environmentally co MD-30 also safeguards All iM3 Advantage hand

F1 Advantage H

Spray Cap for High Speed and Low

Spray Cap

## iM3 Small Animal Dental Equipment

## 1:1 LED Surgical Straight Nose Cone - Redesigned

• The 1:1 LED Advantage fits any low speed motor or electric micromotor • Self generated LED bright white daylight, produced by the drive shaft of the low

Large illumination area providing outstanding visibility

Coolant delivery tube, removable/adjustable sterile coolant to any bur length



"The light on the 1:1 LED handpiece is magic. I used it last week on a few wombats from Australia Zoo (we see them every 3 months) and it made our job a lot easier.

Previously as you know we tried to use things like overhead lights, otoscopes and even the light source from an endoscope. Nothing compares to this device"

Dr Aaron Forsayeth BVSc (Hons) MACVSc (Veterinary Dentistry) Australian Veterinary Dental Society President

Straight Nose Cone 1:1 ratio (blue) uses HP Burs for rabbit and

Ideal for grinding molar teeth in rabbits using long HP diamond burs. Fits any Low Speed air motor or electric micromotor.

Handpiece Oil – 30 ml	Cat. No MD30
rt maintenance intervals = longer lifespan tic oil specifically designed for the increased v dentistry xic constituents is physiologically harmless compatible s the durability of plastic parts, coatings and seals etc. adpieces must be lubricated with MD-30	
Handpiece Oil - Aerosol	MD400
v Speed Handpieces	MD4010

for 1:1 SNC LED and Low Speed Attachments


Cat. No 16860

MD4020

Cat. No

L7800

The Ultrasonic difference - Amplitude, Frequency and Tip Motion.

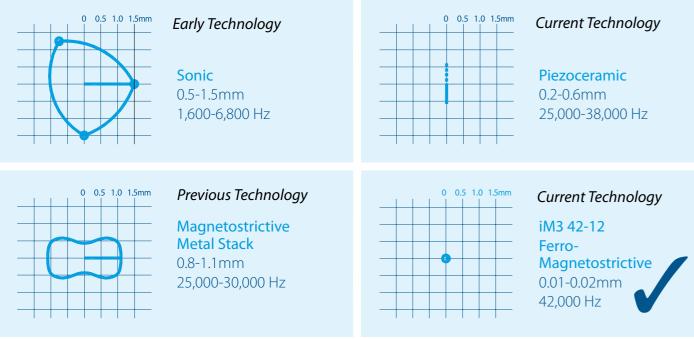

Smaller tip motion (movement) means less enamel damage.

#### 42-12 Ultrasonic scaler (42,000 cycles) iM3

The 42-12 scaler is unique as it uses the high ultrasonic frequency to remove the plaque and calculus not the physical scraping of the tooth as do the scalers mentioned above. The 42-12 tip is totally round and rotates in a circular motion 0.02mm causing little or no enamel damage. The tip has no sharp sides like the Piezo or Cavitron style scaler. The 42-12 scaling tip is made from titanium and is fitted with a ferrite rod (not metal). The tip generates little or no heat and can be safely used subgingivally up to 13mm.

The big advantage of the 42-12 is the speed at which a prophylaxis can be preformed and the ability to use safely subgingivally. Many users report that scaling time is reduced by half when compared to previous scalers used.

## Not all ultrasonic scalers are the same!

A survey of 4737 AAHA Members in the USA - of which 667 responded completing a comparison of the 5 major ultrasonic dental scalers in 15 categories - iM3 came in first in every category, scoring a 90% overall satisfaction. TRENDS MAGAZINE

Unlike other ultrasonic tooth scalers that use a "jack hammer" approach to tooth scaling, the 42-12 adopts a more gentle technique with its revolutionary rotational tip movement, and ultrasonic removal of the plaque and calculus.

#### 42-12 Ultrasonic Scaler

The 42-12 rapidly removes heavy calculus deposits from all tooth surfaces and is safe to use in periodontal pockets up to 13 mm deep.

The 42-12 achieves a highly polished finish on the tooth with less enamel damage by using an exceptionally high ultrasonic frequency of 42,000Hz, coupled with an extremely low amplitude of the titanium scaler tip's circular movement of 0.02 mm • Power is permanently set at a high 42,000 Hz frequency (cycles) yet is safe to use in periodontal pockets up to 13 mm deep. No other veterinary ultrasonic scaler can do this safely because they produce too much heat

- coolant mist which is adjustable

- Autoclavable handpiece and tips

#### iM3 42-12 Ferrite Rod and Scaler Tip Insert

Ferrite Rod For the iM3 42-12 Scaler Tip Insert

Universal Tip

Perio Tip

Thinline straight tip is ideal for felines and small mouths.

#### Tip Wear Guide

Use the 42-12 Ultrasonic Scaler Tip Wear Guide to check whether your ultrasonic scaler tip needs to be replaced.

- ultrasonic frequency.



Lightweight handpiece

with finger on/off control

· Lightweight handpiece with finger on/off control right where you need it, eliminating troublesome foot control altogether

The rotational tip movement and high ultrasonic frequency create a fine

• Rapidly removes plaque and calculus from all tooth surfaces, while using a lighter touch on the enamel than required of other ultrasonic scalers • Unit supplied with universal and perio titanium tips and spare ferrite rod.

Titanium tips last 3-4 times longer than stainless steel scaler tips

## Faster, easier to use, quieter; less tooth damage. In a class of its own compared to anything else 42-12 USER UTAH, USA

1. Tip wear will significantly reduce scaler performance as it will alter the

2. The ferrite rod must be tight - a loose or broken rod will cause poor or no vibration. 3. A correctly functioning tip used on a power setting of 3-4 lights will

produce a coolant mist around the tip.

Cat. No 06979

06840

06850

06845

www.im3vet.com

Cat. No 06860





The Ultra LED handpiece is impressive due to its lightweight, ergonomic design and an LED ring made up of 5 light emitting diodes. The easy to dismantle handpiece, including the LED ring is sterilisable.

#### 3 universal and 1 perio tips are available for the Ultra LED scaler. The tip wear guide is provided free of charge and helps identify a worn tip that needs replacing.

Replacement Ultra LED Scaling Tips					
U801U	Universal Tip, suitable for coarse depuration of supragingival calclus				

U802U	Universal Tip, for fine depuration of supravingival calclus
U803U	Universal Tip, to remove supragingival calcus, suitable for interdental cleaning

U801P Perio Tip, removal of subgingival deposits

Piezo ultrasonic scalers have more power and generate less heat than cavitron style scalers. The tip action is linear and the movement of the tip is used to remove plaque and calculus by physically scraping the tooth's surface. The tip should never be used on its point as it will work like a jackhammer (wood pecker) and create a hole in the tooth. The main disadvantage of the Piezo is still the heat generated, enamel damage and the limitations of scaling 1-2mm safely subgingivally.



#### iM3 P6 Piezo Ultrasonic Scaler

Cat. No U2100

The newly designed Piezo P6 dental scaler is compact in size and high in quality. The P6 is designed and manufactured to meet the most demanding scaling environment. System is portable and lightweight weighing only 3.52 lbs; it comes with a built in power setting dial (1-20) that allow users to set power output at their desired level of operation. Both the handpiece and tips are autoclavable and Satelec® compatable.

Technical specifications are as follows:

- · Power Supply: Available in either 115Volts or 230 Volts -
- Working Freq.: 30,000 Hz (cycles) individual stand alone unit.
- · System includes 4 tips and 1 tip wrench.
- Connects directly to the water system of any iM3 air driven dental machine.
- Pump up water tank with fittings available for stand alone use.
- Dimension: 7.68"(L)x3.38"(W)x3.27"(H) Weight: 3.52 lbs.
- Warranted, sold and serviced by iM3.




iM3 P6 LED Piezo Ultrasonic Scaler

Cat. No U2200

The P6 LED Piezo includes all of the great features on the standard P6 as well as a bright white LED ring at the end of the handpiece illuminating the scaling tip and working area.

The LED light has a built in delay keeping the LED light on after you release your foot from the control pedal keeping the working area illuminated. The system includes 4 scaling tips that are autoclavable.

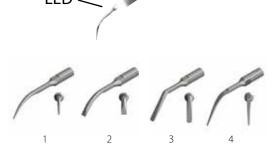

LED

#### **Replacement Piezo Scaling Tips** Suitable for P4 and P6 Piezo (BS)

1. BS-1 Pointed tip (recommended)
2. BS-2 Square Tip
3. BS-3 Round tip
4. BS-P Perio Tip

iM3 Code U2002 iM3 Code U2013 iM3 Code U2014 iM3 Code U2015

Cavitron Style Scaling Tip Insert 25k - Universal

Cat. No C2500

4

## iM3 Ultra LED

The new iM3 Piezo ultrasonic scaler with LED light technology has been designed for use in the areas of scaling, periodontology, endodontics, and micro-preparation. Thanks to the Ultra LED Piezo Scaler, you can be sure to achieve the best results in all applications.

The ergonomic design of the handpiece and the light weight make it possible for you to work for longer without becoming tired. The handpiece itself can be disinfected and sterilized.

#### Your advantages at a glance

- Allows you to work in daylight guality light
- Reflects the animals tooth colour exactly
- Allows you to see contrasts much better
- LED ring with 5 light emitting diodes that remain on for 30 seconds to assist with oral examination
- Ergonomic design
- · Perio and universal scaling tips and holders supplied Inclusive water filter
- Handpiece is sterilizable up to 135 °C
- Tip changer and tip cards included
- 2 year warranty with unlimited iM3 support.


have been using the new Ultra LED piezo ultrasonic scaler for several months now. It is a very powerful and efficient scaler. I find the LED function especially helpful in the hard to see and reach areas of the mouth, the areas such as the back molar areas

Erik Christensen, CVT Arizona Veterinary Dental Specialists



#### 5 LEDs

Light emitting diodes (LEDs) result in robust, sterilisable light that barely heats up and are shock-resistant. LED's last much longer than traditional light sources.

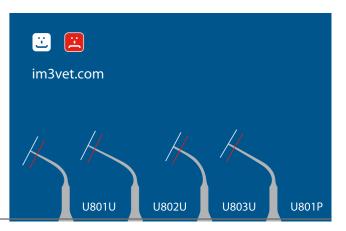

Cat. No 118000

4

## A major advancement in tooth extraction brought to you by iM3!

# **Vet-Tome**



## I love my Vet-Tome!

My extractions are more elegant and less stressful. With the aid of a thin blade the Vet-Tome cuts the periodontal ligament allowing me to get most teeth out with minimal trauma."

> Jean Battig, DVM, DAVDC Animal Dental Clinic NW Advanced Dentistry & Oral Surgery

www.animaldentalclinicnw.com

The Vet-Tome improves the tooth removal process for both vets and their patients by providing greater control during extractions. The Vet-Tome is an automated periotome that offers extremely precise tooth extraction with minimal or no alveolar bone loss.

The surgery is often flapless so the animal experiences reduced pain and swelling. This translates to less time spent extracting teeth and faster recovery time for the animal.

Normal evaluation and closure of the alveolus may still be required after extraction.





# **Vet-Tome**®

The atraumatic extraction system that greatly reduces extraction time and saves bone!

## **VET-TOME BENEFITS**

- Less trauma to the patient
- Speeds up the procedure
- May eliminate or reduce the need • to remove buccal bone
- Allows work in tighter spaces
- Cuts paths for luxators and winged elevators
- Decreases hand fatigue



iM3 is excited to launch the new Vet-Tome, an automated periotome with foot pedal operation that offers precise tooth extraction with minimal or no alveolar bone loss and less trauma. The electric Vet-Tome uses replaceable ultra-thin stainless steel tips/blades coupled with the mechanical in and out action generated by the handpiece to cut the periodontal ligament similar to a luxator (hand instrument). The mechanical action and thin flexible blade allows easy insertion into the periodontal ligament (PDL) and alveolus (tooth socket) coupled with the side to side movement by the operator cutting the PDL.

The Vet-Tome will not work for periodontal ligaments that are mineralized, i.e. it won't go through mineral/bone, but it does an excellent job on the canine teeth (and others) which can be very hard, especially lower cat canine teeth which can lead to jaw fracture. The Vet-Tome is not a replacement for good extraction technique and flap preparation but certainly a must have for any veterinary clinic preforming routine extractions.

V2000 Vet-Tome extraction system Complete with 2 standard and 2 heavy S/S tips.

#### V2001 Vet-Tome standard flat tip - V1 (5 pack) 37mm long

The Vet-Tome standard flat tip is a thin blade designed to follow (bend and curve) the tooth root and cut the periodontal Ligament.

#### V2002 Vet-Tome Heavy tip - V2 (5 pack) 40mm long

The Vet-Tome heavy tip is thicker and slightly longer than the standard tip (V2001) designed to be more robust and therefore last longer. The disadvantage of the heavy tip is it does not bend and curve and should only be used on straight tooth roots.

Cat. No V2000

V2001

V2002

Cat. No

L6000

## LED Curing Light

## New Generation Fiber Type LED Curing Light Efficient, Powerful and Reliable

Our new generation Fibre Optic LED curing light features a lightweight, ergonomic handle, designed to ensure easy access to the oral cavity. Most importantly this unit can be battery (cordless) or electrically operated and able to be conveniently located chair side.

- Cordless curling light with rechargeable base
- Can be used with power sourse connected directly to HP
- Light intensity 1000x/cm
- 450-470 Wavelength
- Weight 70g, 220mm in length
- Depth of cure, 3mm/10 seconds
- Long Use-life LED, 50,000 hours
- Compatble with all bonding systems
- 3 modes: fast, ramp or pulse

5 - 42



Prior to experiencing your handpieces, I had no idea how much easier a surgical extraction could be! iM3 USER, USA

#### iM3 PL3 polisher

The iM3 PL3 polisher is a user-friendly, stylish, and compact system. The polisher features a built-in digital display for viewing motor rotation speed and offers speed control of 0 to 30,000 RPM. The PL3 comes equipped with Automatic Feedback Control and powerful torque.

The iM3 PL3 Includes; Micro-motor, iM3 1:1 straight nose cone, foot pedal and 10 iM3 oscillating disposable prophy cups and paste.

- Dimensions (L x W x H): 19.2cm x 15.5cm X 8.5cm
- Hand piece cable length: 2.5m

- Max Polisher Speed: 30,000 RPM
- iM3 HP dental Burs and contra attachments available
- Warranted, sold and serviced by iM3.

iM3 SP6 P6 Piezo Scaler and Micromotor

The iM3-SP6 combines the latest Piezo P6 technology and micromotor in a compact user friendly design. The power setting for both the scaler and micromotor are easily adjusted via the digital LCD display and the micromotor rotates up to 30,000 rpm. The iM3-SP6 is supplied with a 1:1straight nose cone allowing the system to be used for polishing with the prophy angles supplied and the speed turned down to approximately 4,000 rpm.

 ${\rm HP}$  burs available from iM3 can be used in the straight nose cone for cutting and sectioning teeth.

#### iM3-SP6 is supplied with:

- 2 universal scaling tips
- Footpedal control on/off
- Micromotor with E fitting
- iM3 1:1 straight nose cone (Blue) starter polishing pack including prophy angles and paste, water line with male QD connector.



Cat. No U2010

Cat. No

L4200

# www.im3vet.com

tor Cat. No L4300

Optional Pump up water tank with fittings for standalone use is available from iM3.





# M3

The GS Deluxe dental machine has a two bottle coolant system. The clear bottle is filled with distilled water, while the brown bottle protects the UV light-sensative CLS solution. CLS will help minimise cross contamination when used between patients. The spearmint smell will mask mouth odours.

> 3 Way syringe - ideal for washing away prophy paste and debris



expect from iM3.

iM3 Advantage 2pc low speed polisher

3 way air water

the oral cavity

syringe for rinsing

overheat

5 - 7



handpiece for fast sectioning

#### iM3 GS Deluxe Dental Unit

The GS Deluxe is a complete dental system. It includes a high quality Advantage push button high speed handpiece (400,000 rpm) with triple coolant spray for fast sectioning of teeth, a 3-way air water syringe with the ability to select either water, CLS, air or a mist and a low speed handpiece for polishing using iM3 R&R oscillating prophy angles that won't wind up muzzle hair.

The three way (air, water and mist) syringe is a useful instrument for washing away debris from the mouth. The tip is removable and autoclavable and can be rotated 360 degrees. The use of CLS in the dental machine allows for easy application of chlorhexidine into the mouth during dental procedures. Stainless steel instrument tray with silicon mat provides added storage and the ideal location for your iM3 dental instruments.

The iM3 42-12 ultrasonic scaler (pictured p35) or any brand of scaler can be connected to the water and power supply of the GS Deluxe dental machine. The dental machine is self-contained and only one power cord is required to power the machine and ultrasonic scaler.

iM3 GS Deluxe (non LED) with oil free compressor

Cat. No U3012

www.im3vet.com

5 - 45

LED Advantage Swivel High Speed


Handpiece. Excellent ergonomics.



iM3 GS Deluxe LED Dental Unit

The iM3 GS Deluxe LED sets the standard for veterinary dental machines. The GS Deluxe LED is a complete dental system with advanced features not found on any other dental machine.

Quality handpieces coupled with a 5 year dental machine and oil free crompressor warranty provides peace of mind.

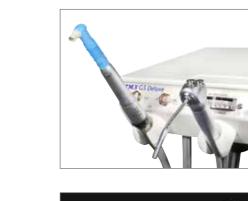

TM 3

Features and Benefits of the GS Deluxe LED Dental Unit The GS Deluxe LED is a complete dental system with advanced features not found on any other dental machine. The industry first self generating LED light Advantage high speed swivel handpiece removes the need for troublesome fibre optic wiring, expensive krypton bulbs, switches and transformer. The low speed handpiece and nose cone swivels for comfort and reduces strain on the hand. The GS Deluxe is fitted with two coolant bottles. One for the CLS solution and the other for water. Maintenance instructions are printed on the back of the machine right where you need them.

handpiece

The SW machine is supplied with a Low Speed Motor (with forward and reverse) and 1:1 LED contra angle handpiece with water. Recommended by some European specialists for removing bone.

The SW machine also includes the 4:1 straight nose cone for polishing (L6250 p32).

Stainless steel instrument tray provides added storage and the ideal location. The stainless steel instrument tray is height adjustable and removable.

The iM3 designed magnetic bur holder prevents burs from falling on the ground, they are held in place right where you need them. The Corian<sup>®</sup> surface is impervious to bacteria and easy to clean unlike fibreglass and plastics.

revenue.

are required.


iM3 GS Deluxe LED SW (Slow speed with water cooling)

\* Not available in all countries.

Oscillating Prophy angles won't catch muzzle hair. Air water syringe delivering water / CLS solution and air.

Stainless steel adjustable height

instrument tray

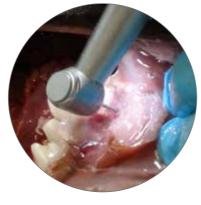

#### **LED High Speed Handpiece Illumination** Light in the mouth is a must!

The new LED Advantage high speed swivel handpiece is the most advanced in the world. It requires no electricity, just air to produce the brightest daylight quality LED white light in the industry with the largest illumination area. Triple coolant spray provides necessary cooling of the tooth and bur. The patented hygiene head helps to eliminate bacterial suckback and cross contamination.

iM3 GS Deluxe LED with Oil Free Compressor

Cat. No U3019

5 - 46

www.im3vet.com

The GS Deluxe LED is supplied with the following high quality handpieces: • High speed swivel LED handpiece • 3 way syringe (air, water and mist) • Low speed

#### NEW GS Deluxe LED SW (Slow Speed Water) Range

In addition to the standard features found on the GSD LED the new SW range includes water on the Low Speed handpiece.

The machine features an adjustable coolant flow valve allowing the user to regulate the amount of coolant from the high speed handpiece.

All iM3 toggles, buttons and auto HP valves can be changed in seconds due to the cartridge valve system used by iM3. This advanced feature allows for the quick cleaning and replacement of valves reducing machine down time and loss of

iM3 developed the first oil free compressor fitted to a dental machine with suction. Oil free technology eliminates the risk of oil contaminates entering the air supply of the dental machine. It also greatly reduces the maintenance regime. The compressor footprint is 51 x 53cm (20 x 21"). The oil free motor generates less heat and is therefore more reliable than conventional oil compressors.

The covered compressor makes cleaning a breeze. iM3's unique technology means the compressor air is cooled on the way to the dental machine. The side of the compressor is left open to allow adequate ventilation, no continuous running fans

Large diameter hospital grade casters make moving the dental machine simple and easy.


Three way (air, water and mist) syringe is a useful instrument for washing away debris from the mouth. The tip is removable and autoclavable and can be rotated 360 degrees. The use of CLS in the dental machine allows easy application of chlorhexidine into the mouth during dental procedures.

#### iM3 GS Deluxe LED with oil free compressor

Cat. No U3019

Cat. No U3019SW www.im3vet.com

| 5 - 47

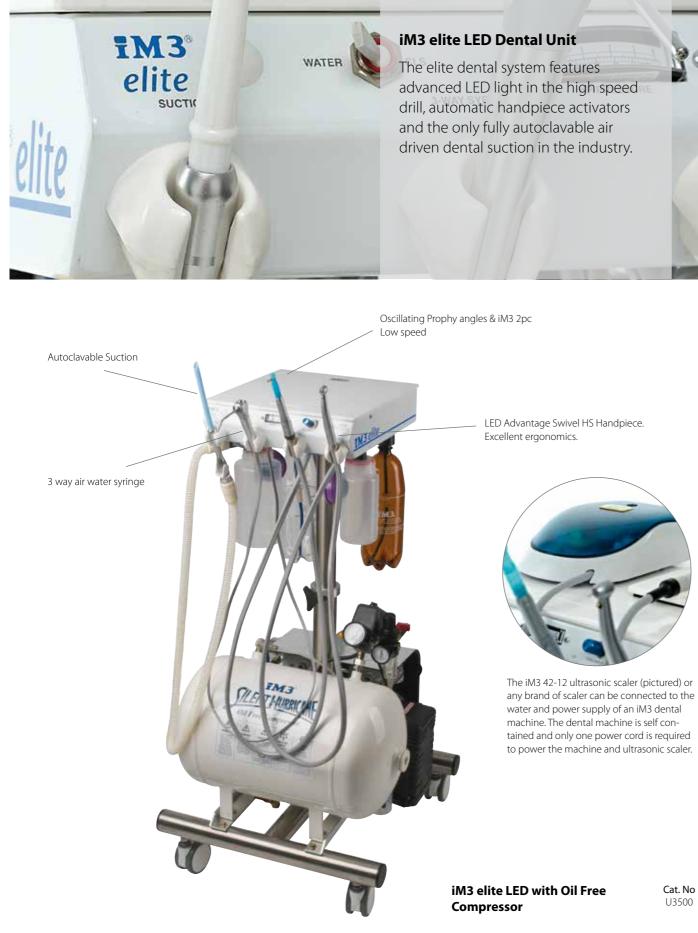



iM3 has developed a variable suction dental unit that is totally autoclavable. The suction tips, handpiece and smooth bore suction tubing can be easily removed and autoclaved.

The suction waste is collected in a 1L container with a replacable filter. Easily accessible and not hidden out of sight! A separate footpedal operates the suction, allowing it's use in conjunction with the high speed drill or 3 way water/air syringe. Because it's fully autoclavable it can also be used in surgery.

The quality of the handpieces and iM3's unique features ensures the elite dental machine is regarded as the number one in the industry. The iM3 elite comes with: Suction • 3 way syringe (air, water and mist) • Low Speed handpiece • High speed LED handpiece with swivel - excellent ergonomics (see page 30)

COOLANT FLOW LS

handpiece.

www.im3vet.com

5 - 48

ntal Units

#### Features and Benefits of the iM3 elite Dental Unit

The iM3 designed magnetic bur holder prevents burs from falling on the ground. The burs are held in place right where you need them. The Corian® surface is impervious to bacteria and easy to clean unlike fibreglass and plastics.

The dental machine chassis is made from aluminium with hospital grade powdercoating. The machine control box and components are covered by a full 5 year warranty. The elite with its built in IEC powerboard and coolant outlet allows for the easy integration of an ultrasonic scaler.

The elite can be supplied with a height adjustable (50cm) long reach (150cm) pneumatic wall arm. The arm allows the machine to be stored against the wall and extended to the end of the wet table for use.

The air line and electricity for the dental machine are contained in the arm, removing unnecessary cables from the ground.

All iM3 machines are fitted with a low speed motor and swivel 4:1 straight nose cone (SNC) for polishing. An array of low speed attachments for rabbit dentistry, orthopaedics and restorative dentistry are available.

Not all automatic handpiece switches are the same! iM3 has developed a pneumatic switch with no electric parts. The switch can be replaced in seconds offering peace of mind. The switch removes the need to manually select the high or low speed

#### iM3 elite LED with oil free compressor

Cat. No U3500 www.im3vet.com


| 5 - 49



#### iM3 Pro 2000 LED Dental Unit

The iM3 Pro 2000 combines all the features iM3 have to offer in a compact self contained dental system. The added storage under the Corian<sup>®</sup> top and the option of having two LED high speed handpieces make this system a must for the busy veterinary practice.

Magnetic bur holder to prevent burs falling on the ground.

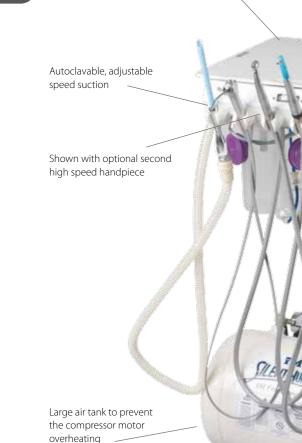

Large hospital grade castors

The Pro 2000 is supplied with an LED Advantage high speed 360 degree swivel handpiece. The swivel allows greater control of the handpiece

Corian® Top, impervious to bacteria, easy to

clean and maintain. Storage underneath

Adjustable height stainless steel stand

and reduces hand fatigue. See P31

Quiet, low maintenance oil free compressor with CE approval

iM3 Pro 2000 LED with Oil **Free Compressor** 

Cat. No U3080

Our technicians fought over who got to use the iM3 so we had to buy a second one!

PRO2000 User, Oxford, OH, USA



Shown with optional second high speed hand piece



handpiece.

floor.



COOLANT

BOTTLE

PRESSURE

COOLANT

OUTLET

www.im3vet.com

5 - 50

ental Units

#### Features and Benefits of the iM3 Pro 2000 Dental Unit

The Pro 2000 is fitted with an LED Advantage high speed handpiece that features 360 degree swivel reducing unnecessary hand fatigue and allows greater control of the

The Pro 2000 comes with the following quality handpieces:

Suction • 3 way syringe (air, water and mist) • Low Speed handpiece

• High speed Swivel LED handpiece • LED light wand 4:1 straight nose cone for polishing and a 1:1straight nose cone for restorative dental attachments, rabbit and orthopaedic surgery using HP burs. All iM3 low and high speed handpieces are autoclavable.

Furthermore the Pro 2000 dental unit is supplied with an additional handpiece cable which gives you the option of having two LED high speed handpieces or an additional low speed handpiece with water cooling.

This high end dental unit is recommended for the busy veterinary practice and it gives the practice all the features iM3 have to offer in a compact self contained dental system.

The Pro 2000 offers a large useable storage space. The smart hinge system means the corian top stays flat when its top is opened. Unlike other dental machines with tops that open up from a hinge at the back, the iM3 hinge system means the items on top do not fall on the

The built in magnetic bur holder prevents the burs from falling on the floor.

The Pro 2000 includes an LED light wand for oral illumination. The light allows you to explore the mouth, see cracks in teeth, while providing additional light in the mouth.

The Pro 2000 is fitted with a two bottle water system with auto pressure release. The advantage of the two bottle system is that it allows the user to choose between CLS (chlorhexidine based solution) and water depending on the procedure.

It's difficult and expensive to send a large (all in one) dental machine away for servicing. The Pro 2000 is built to last. The cartridge and water valves designed by iM3 will keep working well beyond the life of an ordinary valve. Our systems are unique and our engineering is so functional, you will be impressed by the minimal amount of maintenance required to keep your system in top working order.

#### iM3 Pro 2000 LED with oil free compressor

**Cat. No** U3080

5 - 51



5 - 52





#### Features and Benefits of the iM3 Dental Units

The iM3 features make the difference. Our attention to detail and hands on experience in veterinary dentistry means our machines are designed to meet the needs of today's veterinarian. From the general practioner to the advanced veterinary dental specialist, iM3 equipment is the machine of choice for general practice veterinarians, universities, zoos and teaching hospitals throughout the world.

All iM3 air driven dental machines use the same high quality components developed by iM3 and are backed by a 5 year warranty on the dental machine and 1 year on the handpiece.


Pressurised coolant (water or CLS) outlet for any brand of ultrasonic scaler. All maintenance instructions are printed on the back of the dental machine, right where you need them.

All iM3 toggles, buttons and auto handpiece valves can be changed in seconds due to the cartridge valve system used by iM3. This advanced feature allows for the quick cleaning and replacement of valves reducing machine down time and loss of revenue.

Water filters (65 micron) minimise contaminents entering the dental system.

iM3 developed the first oil free compressor fitted to a dental machine with suction. Oil free technology eliminates the risk of oil contaminates entering the air supply of the dental machine. It also greatly reduces the maintenance regime. The compressor footprint is 51 x 53cm (20 x 21") the oil free motor generates less heat and is therefore more reliable than conventional oil compressors. Large diameter hospital grade casters make moving the dental machine simple and easy.

All internal air and water lines are colour coded for easy repair should it become necessary, most of which you can carry out yourself.

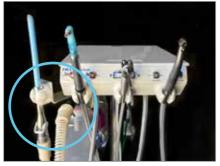

Autoclavable/sterile suction can be retro fitted to existing iM3 dental machines. **Sterile Yankauer Suction sets.** 



iM3 Code S2600





\* All iM3 dental machines can be wall or floor stand mounted with an external air source, using iM3's quiet oil free compressor.

#### Features and Options

It's the iM3 features that make all the difference.

	GS Deluxe Start	GS Deluxe	GS Deluxe LED	elite LED	Pro 2000 LED	SW Range Includes
Standard features on all dental control tops			LED	LED	LED	
Cabling for high and low speed handpieces	Yes	Yes	Yes	Yes	Yes	
Water supply circuit w/ filter, regulator, auto pressure release	Yes	Yes	Yes	Yes	Yes	
2 Bottle water system – CLS and water selector	1 Bottle	Yes	Yes	Yes	Yes	
3-Way air/water syringe with removable autoclavable tip	Yes	Yes	Yes	Yes	Yes	
Quick reference maintenance instructions printed on the machine	Yes	Yes	Yes	Yes	Yes	
Pressurised water outlet for ultrasonic scaler	Yes	Yes	Yes	Yes	Yes	
Suction unit – autoclavable variable suction and bacterial filter, supplied with 5 small and 5 large suction tips	Optional	Optional	Optional	Yes	Yes	
LED light wand					Yes	
Low speed cabling with internal water supply				Optional	Yes	Yes
Second high speed or low speed cabling with coolant flow control					Yes	
4 Socket IEC power strip with circuit breaker	Yes	Yes	Yes	Yes	Yes	
Spare IEC male / female power cord for ultrasonic scaler	Yes	Yes	Yes	Yes	Yes	
Storage compartment	0.1.1					
Stainless steel instrument tray	Optional	Yes	Yes	Optional		
Corian <sup>®</sup> control top	Yes	Yes	Yes	Yes	Yes	
Corian <sup>®</sup> control top with built-in magnetic bur holder	Optional	Yes	Yes	Yes	Yes	
Collection bottle with bacterial filter for handpiece exhaust lines	Optional	Yes	Yes	Yes	Yes	
Colour coded delivery system for ease of service	Yes	Yes	Yes	Yes	Yes	
Cartridge valves to minimise downtime	Yes	Yes	Yes	Yes	Yes	
Powdercoated aluminium and stainless steel construction	Yes	Yes	Yes	Yes	Yes	
Warranty on dental control tops including components	2 years	5 years	5 years	5 years	5 years	
Handpieces and operation						
Toggle to select high & low speed handpieces	Yes	Yes	Yes			
Automatic handpiece activators				Yes	Yes	
Advantage high speed push button handpiece - L7100 (page 31)	H2030	Yes				
LED Advantage high speed handpiece with LED light and 360° swivel fitting – L7300 (page 30)			Yes	Yes	Yes	
Polisher – iM3 Advantage low speed motor and 4:1 straight nose cone (autoclavable) – L6200 and L6250 (page 32)	Yes	Yes	Yes	Yes	Yes	No
LS Motor with Forward/Reverse and speed control. 1:1 LED contra with internal water 4:1 SNC for polishing. L2010, L2004 & L6250				Optional	Optional	Yes
1:1 Straight nose cone for the low speed motor. Used for contra angles, rabbit/rodent dentals and orthopaedic work – L6860	Optional	Optional	Optional	Optional	Yes	
1:1 LED Advantage straight nose cone with LED white light and coolant tube L7800 (page 33)	Optional	Optional	Optional	Optional	Optional	
Warranty on all iM3 handpieces	6 months	1 year	1 year	1 year	1 year	
Mountings and air compressor						
Adjustable height stainless steel floor stand w/ castors	Yes	Yes	Yes	Yes	Yes	
Silent hurricane oil free compressor & stainless steel floor stand	Yes	Yes	Yes	Yes	Yes	
Wall mounted long reach arm, air operated 150cm reach with 50cm travel up and down.	Optional	Optional	Optional	Optional	Optional	
Warranty on mountings and air compressor	2 years	5 years	5 years	5 years	5 years	
Ultrasonic scaler options						
iM3 Piezo Ultrasonic Scaler with 4 tips	Optional	Optional	Optional	Optional	Optional	
iM3 42-12 Ultrasonic Scaler with 2 titanium tips	Optional	Optional	Optional	Optional	Optional	
iM3 Ultra LED with 2 tips	Optional	Optional	Optional	Optional	Optional	

The SW range includes internal water on the low speed handpiece.

iM3 dental machines are made by iM3! We are happy to customize a dental machine to suit your needs.



#### **INTEGRATED X-RAY**

The Revolution 4DC features the latest technology. The full-swivel ball-joint is revolutionary in every sense. The tube head revolves freely around the sphere, so that practically any position can be achieved. An automatic touch-sensitive lock-on/lock-off device enables effortless repositioning of the tube head.



4 .

#### **ATTRACTIVE**

Available in 4 striling colours, the Pro X is designed and manufactured in Europe. It has a timeless and award-winning design. You and your staff will never get tired of looking at and working with Pro X.

5 - 54



#### **PRO X HD DISPLAY**

Pro X HD Display is available with a long double articulated arm that allows you to view the display while you work giving you full and convenient access to all your dental X-rays. Available with either a single arm or an articulated, balanced double arm with handles.



#### SAFE

Hands free user interface, autoclavable handles, easy to clean smooth surfaces and antibacterial paint help you avoid cross contamination.

Integrated suction and water line disinfection systems available.

The Pro X sets the standard for advanced Veterinary Dental treatment centres.

The iM3 Pro X can be configured to your own specification, you can select from a large range of options to create the ultimate dental treatment centre. No longer will you have to make do with "standard" configurations, iM3 will design and customize a system to meet your needs!









#### iM3 Small Animal Dental Equipment



#### **SEE WELL**

See all surfaces of all teeth while sitting well. Choose working position from 9 through 12 o'clock. The Pro X Daylight LED operating lamo offers shadow free light that is always parallel to your viewing direction.



#### **INDIVIDUAL**

There is a Pro X solution to fit your individual needs. Select up to six instruments and other accessories according to your preferences. The Pro X is electronically controlled allowing you to adjust the function of each instrument individually. The instruments are all balanced so that the handpiece cable and actuator arm carries the weight of the handpiece reducing fatigue and allowing for greater control by the user.


#### FOCUS

Operate all unit and instrument functions hands free and maintain full focus on the oral cavity. You do not need to look away from the work area to grab instruments as they are centrally places over the patient and each time an instrument with spray stops, the work area is automatically dried.

Pro X - Dental Treatment Centre

| 5 - 55

## Size Matters – iM3 Dental X-Ray Canine & Feline Positioning Guide

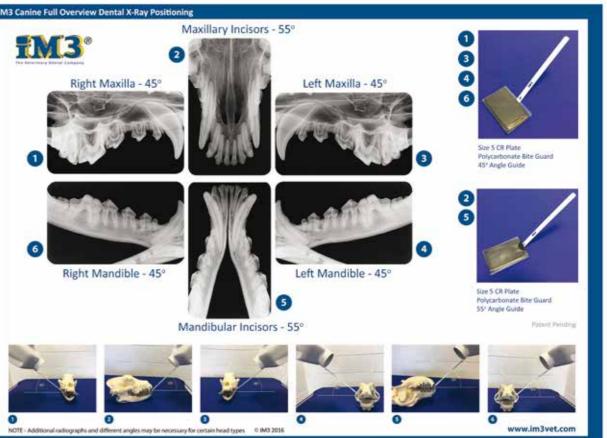

iM3 has designed a new technique to achieve a full mouth series of dental radiographs for every patient. The new technique "Magic 6" allows you to achieve a full mouth series on a medium to large dog or cat in less than 5 minutes with just 6 x-rays. It'is that simple!

With this new technique there are only three angles you need to use, 30°, 45° & 55° and the iM3 Dental X-Ray positioning kit includes these angled guides and plate holders.

Forget trying to use a size 2 DR sensor to x-ray a full mouth in a dog. DR technology has limitations in veterinary dentistry, most notably the small size of the #2 sensor compared to the large size 5 as illustrated. DR sensors are also fragile and expensive to replace. CR plates have a long life and are low cost. The large format CR size 4 and 5 dental plates combined with the positioning guide make it possible to substantially reduce the number of x-rays required for a full mouth series in any animal.



Polycarbonate Bite

Size 4, and Size 5

Image Plates

Protectors for Size 2,



**Canine Positioners** for full mouth X-rays in dogs

## X-Ray Bisecting Angle Guide and Instructions



- Colour coded delivery system aids in fast fault diagnosis or simple replacement by you in your hospital
- Maintenance instructions printed on the chassis of the GS, elite and Pro 2000 machines for guick reference
- New low maintenance oil free quiet compressor designed by iM3 to provide oil free air to the dental machine with handy aluminium cover for storing footpedals when not in use
- Ergonomically designed iM3 low and high speed dental handpieces have been developed specifically for the veterinarian. The design makes changing a prophy cup or dental bur straight forward and will result in a longer working life of the handpiece
- · Manuals available on cd or downloaded via the web www.im3vet.com
- Cartridge valves the toggle valves on all iM3 units use a unique cartridge system allowing the operator to replace the valves quickly, should it become necessary with no down time, another industry first!

Warranties

iM3 products have been manufactured with great care by highly qualified technicians using a wide variety of tests and controls designed to ensure faultless operation.

To maintain a high standard of quality control, every dental system sold by iM3 is thoroughly tested and operated complete with the hand pieces and coolant before it is packaged for delivery. The quality of workmanship and unique design features used to build every iM3 are so good that our range of dental machines now carry a warranty which is setting new standards within the veterinary dental industry. Up to 5 years for our range of air driven dental machines.

At iM3, our reputation and passion is your guarantee of advanced technology, unique features and quality workmanship. We design them, we build them and we service them. We know your iM3 workstation inside and out.

## **Technical Support**

We have had the occasion to call on your

www.im3vet.com

service staff for advice and repairs regarding our dental machine. In every case, your staff have been very knowledgeable and provided assistance with patience and courtesy. You should be congratulated for outstanding service to your clients and we will continue to deal exclusively with you for our dental machine requirements." Peter Roche,

Ferny Grove Veterinary Surgery

#### FREE Technical Support is included with every iM3 Dental Workstation

The relationship with our customers is of the utmost importance to us, to ensure that our products are performing to a standard we're proud to offer. The iM3 Superior Technical Support Team comes FREE with every iM3 dental workstation or iM3 product.


Our technical support team will continue to look after your dental system and all questions regarding handpieces, operational issues or maintenance long after you've taken delivery. We will help you to achieve the most out of your iM3 workstation and make veterinary dentistry a rewarding part of your day.

Call us anytime during business hours on our 'Tech Support Line' and we'll take care of all your needs or email us with your query and we'll get back to you promptly

USA	1 800 664 6348*	info@im3usa.com
AUS	+61 2 9420 5766*	sales@im3vet.com
Ireland	+353 (0)16911277*	sales@im3dental.com
Direct: UI	K +44 (0)1423 224297*	sales@im3dental.com

\*During Normal Business Hours

| 5 - 56

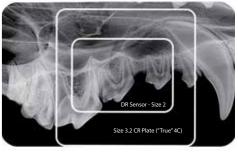

Size 5 CR Plate (Actual Size) - only for the iM3 CR7 VET



Cat. No X8000

The bisecting angle guide provides a simple approach to achieving dental x-ray images. The guide is supplied with notes and diagrams.

www.im3vet.com

6 - 57

X-Ray

"You'd be hard placed to find another piece of equipment that ticks those elusive three boxes of good investment more effectively."



## Dental Radiographs reveal underlying pathology


- Size 2 Feline lower incisors and canines



- - Unique design protects unit and radiographs from dust, hair and particles that can reduce image quality

  - Plate feeder accepts all image plate sizes without the need to change loading cradles - unlike other system Largest range of plate shapes



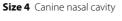

Size 0 Lower mandible of small breed dog

#### The CR7's image quality is unmatched. Having size 4 phosphor plates in addition to smaller plates makes this product superior to any digital computed dental radiography system."



#### Intuitive software for quick and easy operation The synegistically designed and developed CR7 vet software offers a simple, easy to use interface

Using the self-explanatory pictograms, the correct examination parameters can be quickly chosen. Thanks to the flexibility of the image plates, their application is easy and within seconds after the exposure, detailed images can be interpreted onscreen. Full DICOM compliance allows easy export of data to remote applications. DICOM send - optional.

- Unsurpassed high resolution 25 lp/mm, Laser spot 12,5um, 1270dpi
- Scans a range of intra oral plates from 0 5
- Specially designed rabbit image plates available
- Fast processing times (12 seconds size 2)
- Vet dental specific software •
- 10 network licenses
- 2 year warranty
- iM3 factory direct technical support

## It's the range of Plate Sizes and Image Quality that makes the difference!

The CR7 Vet has the largest range of high resolution image plates including sizes 0,1,2,3,4,5 & specially shaped intra oral plates for rabbits!

With DR systems you are often limited to size 2 only! Unlike CR plates that are thin and flexible the size 2 DR sensor is bulky, often too large for cats or small dogs mouths and not large enough for medium to large dogs, resulting in multiple radiographs being required to achieve a full mouth series as well as canine teeth.

## The CR7 is now the No.1choice!

#### CR 7 VET image plate x-ray scanner complete with software

The CR7 Vet system includes:

- CR7 image plate scanner and Vet-Exam software
- 2 year warranty with ongoing iM3 Technical Support



		2	5 CITIX 4 CITI	A/120
currently avail wide with intr	The only scanner	3	2.7 cm x 5.4 cm	X7130
	currently available world-	4	5.7 cm x 7.5 cm	X7140
	wide with intraoral plates for rabbits	5	5.7 cm x 9.4 cm	X7150
	IOLIADDILS	6	Small Rabbit Plate Set	X7165

I have used a number of DR systems in the past, both in veterinary and human practice (Schick, Sirona, Kodak and Genoray), but I would have to say that the results and image quality that I am getting with the iM3 CR7 Vet is the best so far.

The advantages of the CR7 Vet over other DR systems when used in the veterinary environment include a unique range of plate sizes from size 0 up to size 5, which covers all pets from small to large. There is even an intraoral plate for rabbits."

#### **Dr. Anthony Caiafa** BVSc BDSc MACVSc (SA Surgery and Veterinary Dentistry)

References: 2. Verstraete FJ et al. Diagnostic value of full-mouth radiography in cats. Am J Vet Res 1998; 59(6): 692-5.3. Verstraete FJ et al. Diagnostic value of full-mouth radiography in dogs. Am J Vet Res 1998; 59(6): 686-91. †

| 6 - 58

10 Network licences and PACS functionality

· Protective dust cover



**CR7** Vet

Cat. No

X7001

**Dental X-Ray** 

	FIALE NO.	JIZE	IND CODE
	0	2 cm x 3 cm	X7100
	1	2 cm x 4 cm	X7110
	2	3 cm x 4 cm	X7120
The only scanner	3	2.7 cm x 5.4 cm	X7130
currently available world-	4	5.7 cm x 7.5 cm	X7140
wide with intraoral plates for rabbits	5	5.7 cm x 9.4 cm	X7150
or raddits	6	Small Rabbit Plate Set	X7165





Rabbit molars

Intra Oral



Delivers solid profit

Does not need a specialist level of training to deploy

Dr Dave Nicol, MRCVS.

Intraoral radiographs have revealed underlying pathology in: <sup>2,3</sup>

27.8% of dogs where no abnormal findings were found on oral examination 41.7% of cats where no abnormal findings were found on oral examination


#### CR7 Vet. The benchmark in digital X-ray Imaging for veterinary dentistry

- Excellent quality, superior resolution images via Vet-Exam software
- Can be networked acros 10 PCs in the practice
- Full DICOM compliance ensures easy export of data
- Incredibly compact at 23cm x 24cm. Only 6kg.
- Easily removed plate feeder for cleaning
- Precision German manufacture with full 2 year warranty
- iM3 techincal support

Dr. Alexander M. Reiter Dipl. Tzt., Dr med.vet., Dipl. AVDC, EVDC Assoc Professor of Dentistry and Oral Surgery University of Pennsylvania School of Veterinary Medicine

The practice friendly vet specific software allows easy identification of owner and animal specifics, evaluation of image, diagnosis and patient reports complete with radiographs.

6 - 29

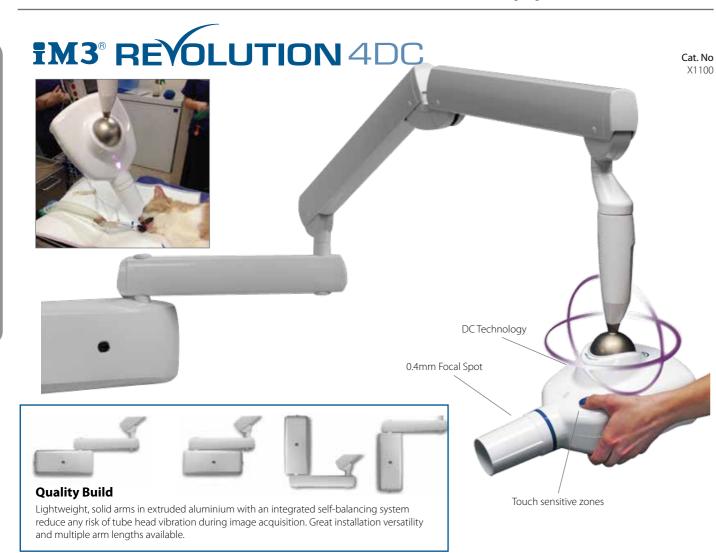

The Revolution 4DC features the latest technology. The full-swivel ball-joint is revolutionary in every sense. The tube head revolves freely around the sphere, so that practically any position can be achieved, including vertical. An automatic touch-sensitive lock-on/ lock-off device enables effortless repositioning of the tube head between exposures.

#### **Comfortable grip**

| 6 - 60

Ergonomically shaped areas on either side of the tube head make rapid positioning easy and effortless to achieve. Right at your finger tips.

#### Touch-sensitive zones

Two conveniently placed touch sensitive zones unlock the tube head for quick and easy positioning. As soon as the operator lets go, the tube head locks into place with great precision.

#### Lightweight and compact

Revolution 4DC is equipped with a 12" X-ray collimator almost entirely integrated into the tube head body. This makes the tube head itself a more compact lightweight component than would otherwise be possible. The unique tube head design, together with the stable aluminium arms, make this unit effortless to manoeuvre.

#### Simple wireless control

One handheld wireless digital control device allows you to program the unit from wherever you are in the surgery. There are no other wall-mounted or wired control panels to worry about. Easy-to-use and easy-to-handle, the control device offers a full range of intuitive exposure scenarios designed to make correct X-ray acquisition a straightforward process. No more complex programming or button-crowded panels. The Revolution 4DC automatically defines the correct exposure by selecting the region of interest.

#### **Precision Imaging**

Cat. No

X1200

Revolution 4DC floor

stand - Optional

Specifically conceived for digital X-Ray imaging, the constant potential tube head combined with the smallest focal spot available for intraoral imaging (0.4mm) guarantees consistently sharp images.

#### **Revolution 4DC Multi-Mode**

Totally automatic choice of the most appropriate combination of technical factors. KV, mA and exposure time are determined according to patient size and anatomical region of interest.


#### **Health and Safety**

Revolution 4DC, with its constant potential, high frequency X-ray generator, significantly reduces radiation dosage compared to traditional X-ray units. Not only is the most harmful, low energy radiation almost completely eliminated thanks to high-efficiency generator, but the embedded 12" collimator minimizes the exposed bodily surface and increases the X-ray parallelism. This leads to better image guality and higher regard for patient's health.

#### **Technical data**

Generator	Constant potential, microprocessor controlled
Operating frequency	145 - 230 KHz (175 KHz typical)
Focal spot	0.4 mm, (IEC 336)
Anode current	6, 7 mA
Voltage (kV setting)	60, 63, 65 kV
Exposure Time	0.010 – 1.000 sec, R10 and R20 scale
Irradiated field	Ø 60 mm round
Total filtration	2.5 mm Al



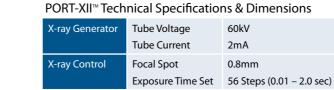

Veight

Dimensions

Main Body

Main Body

2.35kg

197 x 147 x 145 (mm)





#### iM3 Small Animal Dental Equipment



50/60 Hz, 115V ± 10% or 230V ± 15%
self regulating, continuous operation up to
1s/60s
Automatic touch-sensitive lock-on/lock-off
(HyperSphere)
Available in 3 lengths: 40 cm – 60 cm – 90 cm
230 cm, from wa <b>ll</b>
CE 0051, cCSAus, FDA approved

#### Portable Dental X-Ray

Cat. No X1000

The PORT-X II provides total portability combined with the latest digital technology.

• Light and small, enabling total portability and use in any location.

• Small dose of radiation, ensuring patient and operator safety.

• High performance Tube and HFG, providing sharp and clear images.

• Designed for use with both digital sensors and CR image plates.

· Supplied with remote operating switch.


• Fixed kV (60) and mA (2) means the only variable is the exposure time.

• X-ray tube manufactured by Toshiba.

· Complies with National Safety Standards and guidelines.

X-Ray

#### **iM3 Practical Products**

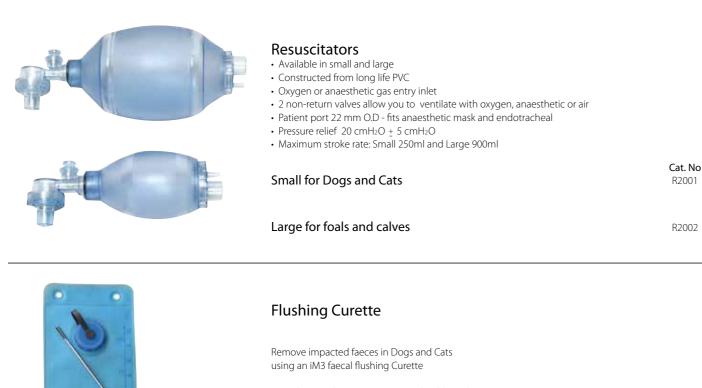

 Stainless Steel construction, autoclavable and a proven design in practice use for over 30 years

- Reduces anaesthesia time
- 2 litre refillable fluid bag
- Gravity flow with a control valve • Allows the optimum solution to be used
- Makes an unpleasant task easy

• Depth markings on curette

T3000 -The complete system

#### **Corporate HQ - Sydney, Australia**



#### iM3 Inc, - Vancouver WA, USA

Cat. No

T3000



#### iM3 Dental Limited - Co Meath, Ireland





#### iM3 Anaesthetic Masks

Unbreakable, see-through low dead space anaesthetic masks. Designed in Australia by iM3 available in 5 sizes with soft diaphragms.

	Cat. No
Anaesthetic Mask - Bird/Rodent	D1940
Anaesthetic Mask - Cat	D1930
Anaesthetic Mask - Small Dog	D1920
Anaesthetic Mask - Large Dog	D1890
Anaesthetic Mask - Long Nose Dog	D1870
5 piece set	D1910

# **DISCOVER THE iM3 DIFFERENCE!**

iM3 is a family owned veterinary Dental Company established in Sydney Australia over 28 years ago.

In that time, iM3 has grown to become a leading veterinary dental company setting the standard in product innovation, industry firsts and continuing education.

From its three locations in Sydney Australia, Vancouver WA, USA and Co Meath, Ireland (Europe), iM3 supplies veterinarians throughout the world with high quality yet affordable veterinary dental products.

iM3 continues to design, develop and manufacture the latest veterinary dental machines, instruments, X-rays and dental consumables. Our large involvement in CE throughout the world means the iM3 staff responsible for the development of new products also know how to use them!

iM3 prides itself on providing the best possible after sales and customer technical support in the industry.

Our dental only focus and the fact that iM3 manufacturers its own dental machines at our three locations means you are dealing with knowledgeable / experienced people who can offer you sound advice.

You can have confidence knowing that leading dental specialists worldwide use and recommend iM3 products.

The iM3 name is your guarantee of guality veterinary specific dental products coupled with the very best aftersales support and warranty in the industry.

Discover the difference!

Yours Sincerely,

#### The iM3 Team!

iM3 is proud to support the ESAVS







FURTHER INFORMATION & REGISTRATION: Phone: +49 6782 23 29 • Fax: +49 6782 43 14 info@esavs.org • www.esavs.org





#### iM3 Pty Ltd - Australia/Asia

21 Chaplin Drive Lane Cove NSW 2066 Australia

Ph: +61 (0)2 9420 5766 Fax: +61 (0)2 9420 5677 Email: sales@im3vet.com www.im3vet.com

#### iM3 Inc. - USA/Canada

12414 NE 95th St, Vancouver, WA 98682 USA

Ph: +1 800 664 6348 Fax: +1 360 254 2940 Email: info@im3usa.com www.im3vet.com

#### iM3 Dental Ltd - Europe

Unit 9, Block 4, City North Business Park, Stamullen, Co. Meath, Ireland

Ph: +353 (0)16911277 Direct: UK +44 (0)1423 224297 Email: sales@im3dental.com www.im3vet.com

iM3 products can be found throughout the world. Please visit the iM3 website for a full list of international distributors.



YOUR AUTHORISED DISTRIBUTOR

Atoomweg 476, 3542 AB Utrecht +31 30 241 42 46 info@ud-vet.nl | info@ud-vet.be www.ud-vet.nl | www.ud-vet.be



All prices exclude Tax. Prices subject to change without notice.